



INSTRUMENT PROCESSING SHEET

Agency Kennedy Space CenterS/N 80-006759Florida Department of
Law EnforcementDate In 04/03/2025 DI Completion Date 4/21/25☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>ALL</u> Date <u>04/03/2025</u> ☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: instrument arrived with a piece of the dry gas regulator hose broken off inside the calibration inlet port <u>04/03/2025 ALL</u>		Quality Checks By <u>SLH</u> Date <u>4/8/2025</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>134</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 103</u> 32 mm <u>0.160</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.246</u> (.228 - .278) 103 mm <u>0.519</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP6292</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP6293</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG325603 09/13/2025</td></tr></tbody></table>		Simulator	Serial #	Lot #/Exp	0.050	MP6291	202406K 06/19/2026	0.080	MP6292	202406L 06/19/2026	0.200	MP6293	202406N 06/20/2026	0.080 DGS	N/A	AG325603 09/13/2025	Flow Calibration By _____ Date _____ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By <u>SLH</u> Date <u>4/8/2025</u> ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other _____ <u>Pulled out broken piece of dry</u> <u>gas hose connector from the</u> <u>calibration inlet port and replaced</u> <u>with a new connector.SLH 4/8/25</u>																																												
Simulator	Serial #	Lot #/Exp																																																													
0.050	MP6291	202406K 06/19/2026																																																													
0.080	MP6292	202406L 06/19/2026																																																													
0.200	MP6293	202406N 06/20/2026																																																													
0.080 DGS	N/A	AG325603 09/13/2025																																																													
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table> Notes/Suggested Service: <u>Difficulty attaching controls to the instrument ports.</u> <u>Repeated 0.200 g/210L stability - see data for note.</u> <u>SLH 4/8/25</u>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection By <u>SP</u> Barometric Pressure ID# <u>28662</u> Gauge <u>1018</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2024-B</u> <table border="1" style="width:100%; border-collapse: collapse;"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table> Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 X2 ☐ Other _____ ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.04.25 12:53:54 -0400 Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2025.05.01 15:16:01 -0400 Tech Review / Date _____ Admin Review / Date _____		Simulator	Serial Number	0.000	MP6289	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293
Simulator	Serial #	Lot #	Expiration																																																												
0.000		N/A	N/A																																																												
0.040																																																															
0.100																																																															
0.200																																																															
0.300																																																															
0.080 DGS	N/A																																																														
Simulator	Serial #	Lot #	Expiration																																																												
0.050																																																															
0.080																																																															
0.200																																																															
0.080 DGS	N/A																																																														
Simulator	Serial Number																																																														
0.000	MP6289																																																														
Interferent	MP6290																																																														
0.050	MP6291																																																														
0.080	MP6292																																																														
0.200	MP6293																																																														

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: KENNEDY SPACE CENTER

Serial Number: 80-006759

Time of Inspection: 10:37

Date of Inspection: 04/21/2025

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.000 / 0.048	0.078	0.196	0.079
0.000	0.000 / 0.048	0.079	0.196	0.079
0.000	0.000 / 0.048	0.078	0.195	0.079
0.000	0.000 / 0.048	0.078	0.196	0.079
0.000	0.000 / 0.048	0.078	0.196	0.079
0.000	0.000 / 0.048	0.079	0.196	0.079
0.000	0.000 / 0.048	0.079	0.196	0.079
0.000	0.000 / 0.048	0.079	0.196	0.079
0.000	0.000 / 0.048	0.079	0.196	0.079
0.000	0.000 / 0.048	0.079	0.196	0.079

Standard Deviations	0.0000 / 0.0000	0.0005	0.0003	0.0000
---------------------	-----------------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

05: Control Outside Tolerance- SIM NOT SEATED PROPERLY.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

04/21/2025
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006759 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006759</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>KENNEDY SPACE CENTER</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>04/21/2025</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>10:37</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Digitally signed by

Shayla Platt

Date: 2025.04.21

12:47:40 -04'00'

Date

04/21/2025

SHAYLA D PLATT,

Department Inspector

FDLE/ATP Form 69 October 2024

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006759
04/08/2025
Software: 8100.27

- Stabilities

#1

Test	g/210L	Time
Air Blank	0.000	13:56
Control Test	INT*	13:57
Air Blank	0.000	13:57
Control Test	0.198	13:58
Air Blank	0.000	13:59
Control Test	0.199	13:59
Air Blank	0.000	14:00
Control Test Stats		
Average	0.1980	
Std Dev	0.0010	
Rel Std Dev(%)	0.5051	

*Interferent Detect


Operator's Signature


Root cause analysis performed
per TSOP - no user or equipment
error determined. Cleared room
and ran fan for 20 minutes prior
to repeating. Set

KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006759
04/08/2025
Software: 8100.27

- Stabilites

#2

Test	g/210L	Time
Air Blank	0.000	14:23
Control Test	0.197	14:23
Air Blank	0.000	14:24
Control Test	0.198	14:25
Air Blank	0.000	14:25
Control Test	0.198	14:26
Air Blank	0.000	14:26
Control Test Stats		
Average	0.1977	
Std Dev	0.0006	
Rel Std Dev(%)	0.2921	

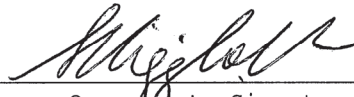

Operator's Signature

KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006759
04/08/2025
Software: 8100.27

- stabilities


Test	g/210L	Time
Air Blank	0.000	14:34
Control Test	0.050	14:35
Air Blank	0.000	14:36
Control Test	0.050	14:36
Air Blank	0.000	14:37
Control Test	0.050	14:38
Air Blank	0.000	14:38
Control Test Stats		
Average	0.0500	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator's Signature

KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000
04/08/2025
Software: 8100.27

-Stabilities

Test	g/210L	Time
Air Blank	0.000	14:29
Control Test	0.079	14:29
Air Blank	0.000	14:30
Control Test	0.080	14:31
Air Blank	0.000	14:31
Control Test	0.080	14:32
Air Blank	0.000	14:32
Control Test Stats		
Average	0.0797	
Std Dev	0.0006	
Rel Std Dev(%)	0.7247	

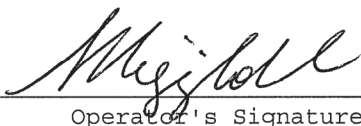

Operator's Signature


KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000
04/08/2025
Software: 8100.27

-stabilities

SN 80-006759

Test	g/210L	Time
Air Blank	0.000	14:41
Control Test	0.081	14:41
Air Blank	0.000	14:42
Control Test	0.080	14:42
Air Blank	0.000	14:43
Control Test	0.080	14:43
Air Blank	0.000	14:43
Control Test Stats		
Average	0.0803	
Std Dev	0.0006	
Rel Std Dev(%)	0.7187	

Dgs

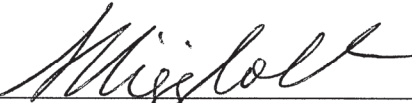

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: KENNEDY SPACE CENTER
Time of Inspection: 13:52

Date of Inspection: 04/08/2025

Serial Number: 80-006759
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08g/210L Test (g/210L) Lot#:_____ Exp:_____	0.20g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:_____ Exp:_____

Number of Simulators Used: _____

Remarks:

BYPASS AI FOR OPERATION COMPLIANCE UNDETERMINED

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

LEANDRA HIGGINBOTHAM

LEANDRA HIGGINBOTHAM

Signature and Printed Name

04/08/2025
Date

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: KENNEDY SPACE CENTER
Time of Inspection: 07:57

Date of Inspection: 04/21/2025

Serial Number: 80-006759
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

Bypassed to operate instrument.

N/A

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Shayla Platt

SHAYLA D PLATT

Signature and Printed Name

04/21/2025
Date