



INSTRUMENT PROCESSING SHEET

Agency FL Highway PatrolS/N 80-007456Florida Department of
Law EnforcementDate In 02/26/2025 DI Completion Date 3/4/2025☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By <u>ALL</u>	Date <u>02/26/2025</u>	Quality Checks	By <u>PN</u>	Date <u>3/4/2025</u>	Flow Calibration	By _____	Date _____																																								
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: <u>drop off, no box</u>			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>223</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 102</u> 32 mm <u>.144</u> (.139 - .169) 36 mm <u>.164</u> (.156 - .190) 53 mm <u>.238</u> (.228 - .278) 103 mm <u>.492</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>30793</u> ☒ Stability Checks			☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																										
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP6292</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP6293</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP6291	202406K 06/19/2026	0.080	MP6292	202406L 06/19/2026	0.200	MP6293	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____																											
Simulator	Serial #	Lot #/Exp																																														
0.050	MP6291	202406K 06/19/2026																																														
0.080	MP6292	202406L 06/19/2026																																														
0.200	MP6293	202406N 06/20/2026																																														
0.080 DGS	N/A	AG429602 10/22/2026																																														
Calibration Adjustment By _____			Department Inspection By <u>PN</u>																																													
Barometric Pressure Gauge _____ ID # _____			Barometric Pressure ID# <u>30793</u> Gauge <u>1013</u> Instrument <u>1013</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2024-B</u>																																													
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6294</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>						Simulator	Serial Number	0.000	MP6294	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293
Simulator	Serial #	Lot #	Expiration																																													
0.000		N/A	N/A																																													
0.040																																																
0.100																																																
0.200																																																
0.300																																																
0.080 DGS	N/A																																															
Simulator	Serial Number																																															
0.000	MP6294																																															
Interferent	MP6290																																															
0.050	MP6291																																															
0.080	MP6292																																															
0.200	MP6293																																															
☐ Post Calibration Adjustment Stability Checks			Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____																																													
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																									
Simulator	Serial #	Lot #	Expiration																																													
0.050																																																
0.080																																																
0.200																																																
0.080 DGS	N/A																																															
Notes/Suggested Service: _____			Digitally signed by Destinee Armstrong Date: 2025.03.05 14:29:10 -05'00'																																													
			Shayla Platt Date: 2025.03.06 11:04:22 -05'00'																																													
			Tech Review / Date _____ Admin Review / Date _____																																													

STABILITY CHECKS

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-107456
03/04/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:52
Control Test	0.048	10:53
Air Blank	0.000	10:53
Control Test	0.048	10:54
Air Blank	0.000	10:54
Control Test	0.048	10:55
Air Blank	0.000	10:56
Control Test Stats		
Average	0.0480	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-107456
03/04/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	10:57
Control Test	0.079	10:57
Air Blank	0.000	10:58
Control Test	0.079	10:59
Air Blank	0.000	10:59
Control Test	0.079	11:00
Air Blank	0.000	11:00
Control Test Stats		
Average	0.0787	
Std Dev	0.0006	
Rel Std Dev(%)	0.7339	

Operator's Signature

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-107456
03/04/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:02
Control Test	0.198	11:03
Air Blank	0.000	11:03
Control Test	0.198	11:04
Air Blank	0.000	11:05
Control Test	0.198	11:05
Air Blank	0.000	11:06
Control Test Stats		
Average	0.1980	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

Operator's Signature

FHP
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-107456
03/04/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:09
Control Test	0.081	11:09
Air Blank	0.000	11:10
Control Test	0.079	11:10
Air Blank	0.000	11:10
Control Test	0.080	11:11
Air Blank	0.000	11:11
Control Test Stats		
Average	0.0800	
Std Dev	0.0010	
Rel Std Dev(%)	1.2500	

Operator's Signature

DGS

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: FHP

Time of Inspection: 14:14

Date of Inspection: 03/04/2025

Serial Number: 80-007456

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2025 <i>2036</i>	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.049	0.078	0.198 <i>PN 3/4/25</i>	0.081
0.000	0.049	0.078	0.197	0.080
0.000	0.048	0.079	0.197	0.081
0.000	0.049	0.078	0.197	0.081
0.000	0.049	0.079	0.197	0.081
0.000	0.049	0.079	0.198	0.081
0.000	0.049	0.078	0.197	0.081
0.000	0.049	0.079	0.197	0.080
0.000	0.049	0.078	0.198	0.081
0.000	0.049	0.079	0.198	0.081

Standard Deviations	0.0003	0.0005	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5


Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

PHIL NICODEMO
Signature and Printed Name

03/04/2025
Date

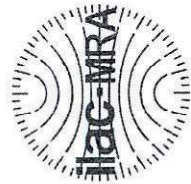

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007456, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007456</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>FHP</u>	0.050 g/210 L	0.004
Calibration Date:	<u>03/04/2025</u>	0.080 g/210 L	0.004
Calibration Time:	<u>14:14</u>	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 October 2024

Issuing Authority: Alcohol Testing Program

03/04/2025

Date

PHIL NICODEMO,

Department Inspector

Service • Integrity • Respect • Quality