



INSTRUMENT PROCESSING SHEET

Agency Hendry CSOS/N 80-006386Florida Department of
Law EnforcementDate In 01/16/2025 DI Completion Date 01/29/2025 Ship P/U H/D CMI EE

Intake	By TDG	Date <u>01/27/2025</u>	Quality Checks	By TDG	Date <u>01/29/2025</u>	Flow Calibration	By	Date																								
<input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE Visual Inspection: <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/ Accessories: <input type="checkbox"/> Power cord <input type="checkbox"/> Printer Cable <input type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: <u>Dropped off. Missing back two feet on DGS shelf.</u>			<input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value <u>187</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column # <u>ATP101</u> 32 mm <u>0.167</u> (.139 - .169) 36 mm <u>0.187</u> (.156 - .190) 53 mm <u>0.257</u> (.228 - .278) 103 mm <u>0.519</u> (.447 - .547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID # <u>26932</u> <input checked="" type="checkbox"/> Stability Checks			Flow Column # _____ <input type="checkbox"/> 5L/min - 17mm <input type="checkbox"/> 15L/min - 53mm <input type="checkbox"/> 30L/min - 103mm <input type="checkbox"/> R-Value _____ <input type="checkbox"/> Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)			Maintenance By _____ Date _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other _____																							
			<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td rowspan="2">MP6286</td> <td>202303K</td> </tr> <tr> <td></td> <td>03/29/2025</td> </tr> <tr> <td>0.080</td> <td rowspan="2">MP6287</td> <td>202303L</td> </tr> <tr> <td></td> <td>03/29/2025</td> </tr> <tr> <td>0.200</td> <td rowspan="2">MP6288</td> <td>202304C</td> </tr> <tr> <td></td> <td>04/05/2025</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG429602</td> </tr> <tr> <td></td> <td></td> <td>10/22/2026</td> </tr> </tbody> </table>			Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K		03/29/2025	0.080	MP6287	202303L		03/29/2025	0.200	MP6288	202304C		04/05/2025	0.080 DGS	N/A	AG429602			10/22/2026			
Simulator	Serial #	Lot #/Exp																														
0.050	MP6286	202303K																														
		03/29/2025																														
0.080	MP6287	202303L																														
		03/29/2025																														
0.200	MP6288	202304C																														
		04/05/2025																														
0.080 DGS	N/A	AG429602																														
		10/22/2026																														

Calibration Adjustment	Department Inspection																																																												
By _____ Barometric Pressure Gauge ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			By <u>TDG</u> Barometric Pressure ID# <u>26932</u> Gauge <u>1021</u> Instrument <u>1021</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2023-B</u> <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP6284</td> </tr> <tr> <td>Interferent</td> <td>MP6285</td> </tr> <tr> <td>0.050</td> <td>MP6286</td> </tr> <tr> <td>0.080</td> <td>MP6287</td> </tr> <tr> <td>0.200</td> <td>MP6288</td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288
Simulator	Serial #	Lot #	Expiration																																																										
0.000		N/A	N/A																																																										
0.040																																																													
0.100																																																													
0.200																																																													
0.300																																																													
0.080 DGS	N/A																																																												
Simulator	Serial #	Lot #	Expiration																																																										
0.050																																																													
0.080																																																													
0.200																																																													
0.080 DGS	N/A																																																												
Simulator	Serial Number																																																												
0.000	MP6284																																																												
Interferent	MP6285																																																												
0.050	MP6286																																																												
0.080	MP6287																																																												
0.200	MP6288																																																												
Notes/Suggested Service: _____ _____ _____ _____ _____	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use <div style="display: flex; justify-content: space-between;"> <div> Destinee Armstrong <small>Digitally signed by Destinee Armstrong Date: 2025.01.31 08:14:08 -05'00'</small> </div> <div style="text-align: center;"> Shayla Platt <small>09:25:07 -05'00'</small> </div> <div> Digitally signed by Shayla Platt Date: 2025.02.03 </div> </div>																																																												
	Tech Review / Date _____ Admin Review / Date _____																																																												

Stability Checks

0.05g/210L 0.047 to 0.053	0.08g/210L 0.077 to 0.083	0.20g/210L 0.194 to 0.206	DGS 0.08g/210L 0.077 to 0.083																																																																																																																																																
<p>HENRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 01/29/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:23</td></tr> <tr><td>Control Test</td><td>0.049</td><td>12:23</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:24</td></tr> <tr><td>Control Test</td><td>0.050</td><td>12:24</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:25</td></tr> <tr><td>Control Test</td><td>0.050</td><td>12:26</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:26</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0497</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1625</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	12:23	Control Test	0.049	12:23	Air Blank	0.000	12:24	Control Test	0.050	12:24	Air Blank	0.000	12:25	Control Test	0.050	12:26	Air Blank	0.000	12:26	Control Test Stats			Average	0.0497		Std Dev	0.0006		Rel Std Dev(%)	1.1625		<p>HENRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 01/29/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:30</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:30</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:31</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:32</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:32</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:33</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:33</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0790</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	12:30	Control Test	0.079	12:30	Air Blank	0.000	12:31	Control Test	0.079	12:32	Air Blank	0.000	12:32	Control Test	0.079	12:33	Air Blank	0.000	12:33	Control Test Stats			Average	0.0790		Std Dev	0.0000		Rel Std Dev(%)	0.0000		<p>HENRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 01/29/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:16</td></tr> <tr><td>Control Test</td><td>0.198</td><td>12:17</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:17</td></tr> <tr><td>Control Test</td><td>0.196</td><td>12:18</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:18</td></tr> <tr><td>Control Test</td><td>0.196</td><td>12:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:20</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1967</td><td></td></tr> <tr><td>Std Dev</td><td>0.0012</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.5871</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	12:16	Control Test	0.198	12:17	Air Blank	0.000	12:17	Control Test	0.196	12:18	Air Blank	0.000	12:18	Control Test	0.196	12:19	Air Blank	0.000	12:20	Control Test Stats			Average	0.1967		Std Dev	0.0012		Rel Std Dev(%)	0.5871		<p>HENRY COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006386 01/29/2025 Software: 8100.27</p> <p style="text-align: right; font-size: 2em;"><i>NS</i></p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:12</td></tr> <tr><td>Control Test</td><td>0.079</td><td>12:13</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:13</td></tr> <tr><td>Control Test</td><td>0.078</td><td>12:13</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:14</td></tr> <tr><td>Control Test</td><td>0.077</td><td>12:14</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:15</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0780</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.2821</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	12:12	Control Test	0.079	12:13	Air Blank	0.000	12:13	Control Test	0.078	12:13	Air Blank	0.000	12:14	Control Test	0.077	12:14	Air Blank	0.000	12:15	Control Test Stats			Average	0.0780		Std Dev	0.0010		Rel Std Dev(%)	1.2821	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:23																																																																																																																																																	
Control Test	0.049	12:23																																																																																																																																																	
Air Blank	0.000	12:24																																																																																																																																																	
Control Test	0.050	12:24																																																																																																																																																	
Air Blank	0.000	12:25																																																																																																																																																	
Control Test	0.050	12:26																																																																																																																																																	
Air Blank	0.000	12:26																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0497																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	1.1625																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:30																																																																																																																																																	
Control Test	0.079	12:30																																																																																																																																																	
Air Blank	0.000	12:31																																																																																																																																																	
Control Test	0.079	12:32																																																																																																																																																	
Air Blank	0.000	12:32																																																																																																																																																	
Control Test	0.079	12:33																																																																																																																																																	
Air Blank	0.000	12:33																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0790																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:16																																																																																																																																																	
Control Test	0.198	12:17																																																																																																																																																	
Air Blank	0.000	12:17																																																																																																																																																	
Control Test	0.196	12:18																																																																																																																																																	
Air Blank	0.000	12:18																																																																																																																																																	
Control Test	0.196	12:19																																																																																																																																																	
Air Blank	0.000	12:20																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1967																																																																																																																																																		
Std Dev	0.0012																																																																																																																																																		
Rel Std Dev(%)	0.5871																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:12																																																																																																																																																	
Control Test	0.079	12:13																																																																																																																																																	
Air Blank	0.000	12:13																																																																																																																																																	
Control Test	0.078	12:13																																																																																																																																																	
Air Blank	0.000	12:14																																																																																																																																																	
Control Test	0.077	12:14																																																																																																																																																	
Air Blank	0.000	12:15																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0780																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	1.2821																																																																																																																																																		
<p style="text-align: right;">Operator's Signature <i>NS</i></p>	<p style="text-align: right;">Operator's Signature <i>NS</i></p>	<p style="text-align: right;">Operator's Signature <i>NS</i></p>	<p style="text-align: right;">Operator's Signature <i>NS</i></p>																																																																																																																																																

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HENDRY COUNTY SO
Time of Inspection: 15:06

Date of Inspection: 01/29/2025

Serial Number: 80-006386
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.050	0.079	0.196	0.078
0.000	0.049	0.079	0.196	0.078
0.000	0.049	0.079	0.197	0.078
0.000	0.050	0.079	0.196	0.077
0.000	0.049	0.079	0.196	0.078
0.000	0.050	0.079	0.196	0.078
0.000	0.050	0.078	0.197	0.077
0.000	0.049	0.079	0.197	0.078
0.000	0.050	0.079	0.197	0.078
0.000	0.050	0.079	0.196	0.078


Standard Deviations	0.0005	0.0003	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies () does not comply () with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

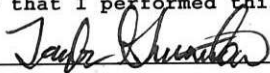

 Signature and Printed Name

TAYLOR D GUTSCHOW

01/29/2025
 Date

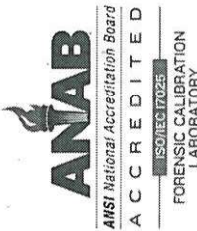

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006386, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006386</u>	UNCERTAINTY* ±	
Owning Agency:	<u>HENDRY COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>01/29/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>15:06</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

01/29/2025

Date TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality