

Agency: KENNEDY SPACE CENTER				Instrument Serial Number: 80-006759					
Date In: 4/10/2026				DI Completion Date: N/A					
				☐ Ship ☐ P/U ☐ H/D ☒ CMI ☐ EE					
Intake By: WKP Date: 4/10/2026				Quality Checks By: SLH Date: 4/30/2026		Flow Adjustment By:			
☒ Annual ☐ Dropped Off ☐ Registration ☐ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☒ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:				☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: 116 ☒ Flow Verification (L/s) Flow Column #: ATP 102 32 mm 0.140 (.139-.169) 36 mm 0.156 (.156-.190) 53 mm 0.230 (.228-.278) 103 mm 0.496 (.447-.547) ☒ Barometric Pressure Check Gauge ID #: 28427 Gauge: 1007 Instrument: 1008 ☒ Stability Checks		Flow Column #: ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: ☐ Post Adjustment Verification (L/S) Flow Column #: 32 mm (.139-.169) 36 mm (.156-.190) 53 mm (.228-.278) 103 mm (.447-.547)			
				Simulator		Serial #		Lot#/Exp	
				0.050					
				0.080					
				0.200		MP6293		202406N	
								06/20/2026	
				0.080 DGS		N/A			
Optical Bench Adjustment By:				Department Inspection By:					
Barometric Pressure Gauge: ID#:				Barometric Pressure ID#:					
Simulator		Serial #		Lot #		Expiration			
0.000				N/A		N/A			
0.040									
0.100									
0.200									
0.300									
0.080 DGS		N/A							
☐ Post Optical Bench Adjustment Stability Checks				Gauge: Instrument:					
Simulator		Serial #		Lot #		Expiration			
0.050									
0.080									
0.200									
0.080 DGS		N/A							
Gauge ID #:				Instrument:					
Gauge:				Instrument:					
Notes/Suggested Service:									
Did not continue with stabilities due to the 0.200 fluctuating results after an ambient fail. In addition, the span noise (in DVM monitoring) was extremely fluctuating on the 9um detector. SLH 4/30/2026 Flow verification was performed on 5/4/26. SLH				☐ Instrument Complies with Chapter 11D-8, FAC					
				☒ Instrument Does Not Comply with Chapter 11D-8, FAC					
				☐ Return to/Place into Evidentiary Use					
				☒ Remain Out of Evidentiary Use					
				☐ Conduct an Agency Inspection Before Evidentiary Use					
Wen-Chi Pierson				Kaitlyn Spearin					
Digitally signed by Wen-Chi Pierson Date: 2026.06.30 08:55:36 -04'00'				Digitally signed by Kaitlyn Spearin Date: 2026.07.15 12:14:13 -04'00'					
Tech Review				Admin Review					

All printed or downloaded copies are uncontrolled
Confirm version prior to use

Document ID: 15246
Version: 2

Issued By: Forensic Quality Manager
Issue Date: December 22, 2025

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: KENNEDY SPACE CENTER

Time of Inspection: 11:11

Date of Inspection: 04/30/2026

Serial Number: 80-006759

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08g/210L Test (g/210L) Lot#:_____ Exp:_____	0.20g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:_____ Exp:_____

Number of Simulators Used: _____

Remarks:

BYPASS AI, COMPLIANCE UNDETERMINED

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



LEANDRA HIGGINBOTHAM

Signature and Printed Name

04/30/2026
Date

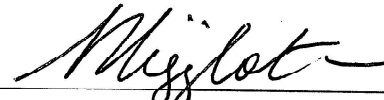
KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006759
04/30/2026
Software: 8100.27

Test	g/210L	Time
Air Blank	AMB*	11:32
Air Blank	AMB*	11:32

*Ambient Fail

0.2009/210L
SIM

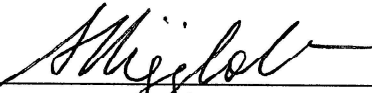
SIM = simulator


Operator's Signature

Ran fan and repeated SUT

KENNEDY SPACE CENTER
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-006759
04/30/2026
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:34
Control Test	0.178	11:35
Air Blank	0.000	11:36
Control Test	0.192	11:36
Air Blank	0.000	11:37
Control Test	0.208	11:37
Air Blank	0.000	11:38
Control Test Stats		
Average	0.1927	
Std Dev	0.0150	
Rel Std Dev(%)	7.7912	


Operator's Signature

Variable results,
Span noise
on the 9mm
detector
fluctuating.
Will send to
CMT.
SWA

Return Material Authorization

Ship to:☒ CMI, Inc.☐ Enforcement ElectronicsShipment to repair facility authorized by: Rodney Snow on 5/28/2026Items Returned: ☒ Instrument ☐ Supplies ☐ Other Describe: _____Instrument Model: Intoxilyzer 8000 Serial Number: 80-006759

Bill To Address:

Kennedy Space Center

Ship to Address:

FDLE Offsite Mail Facility
c/o FDLE HQ
813 B Lake Bradford Rd
Tallahassee, FL 32304

Reason for Return:

Instrument stabilities had ambient fail, variable results for the control, and the span noise displayed large fluctuations. Please also look at the pressure sensor due to value being close to 100.

I require an estimate **BEFORE** any repairs will be authorized and/or conductedPlease contact: Rodney Snow and Terry Quinn Phone #: 321-861-6644 (Snow)Email: Both want to be emailed rodney.f.snow@nasa.gov and terry.a.quinn@nasa.govATP Contact Name: LeAndra Higginbotham ATP Email: leandrahigginbotham@fdle.state.fl.us