



INSTRUMENT PROCESSING SHEET

Agency OCALA PDS/N 80-005454Florida Department of
Law EnforcementDate In 9/2/2025DI Completion Date 09/04/2025☐ Ship ☒ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By <u>WKP</u>	Date <u>9/3/2025</u>	Quality Checks	By <u>SLH</u>	Date <u>09/03/2025</u>	Flow Calibration	By _____	Date _____																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☒ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: <u>Instrument missing two</u> <u>top dry gas shelf feet. WKP</u> <u>9/3/2025</u> <u>Dropped by agency.SLH 9/3/25</u>			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>150</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP102</u> 32 mm <u>0.140</u> (.139 - .169) 36 mm <u>0.156</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28427</u> ☒ Stability Checks			☐ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202406K 06/19/2026	0.080	MP5089	202406L 06/19/2026	0.200	MP5090	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	<table border="1"><thead><tr><th colspan="2">Maintenance By _____ Date _____</th></tr></thead><tbody><tr><td colspan="2">☐ Battery Replacement</td></tr><tr><td colspan="2">☐ Dry Gas Regulator Replacement</td></tr><tr><td colspan="2">☐ Breath Tube Replacement</td></tr><tr><td colspan="2">☐ Other _____</td></tr></tbody></table>			Maintenance By _____ Date _____		☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____				
Simulator	Serial #	Lot #/Exp																																		
0.050	MP5088	202406K 06/19/2026																																		
0.080	MP5089	202406L 06/19/2026																																		
0.200	MP5090	202406N 06/20/2026																																		
0.080 DGS	N/A	AG429602 10/22/2026																																		
Maintenance By _____ Date _____																																				
☐ Battery Replacement																																				
☐ Dry Gas Regulator Replacement																																				
☐ Breath Tube Replacement																																				
☐ Other _____																																				
Calibration Adjustment By _____			Department Inspection By <u>SLH</u>																																	
Barometric Pressure Gauge _____ ID # _____			Barometric Pressure ID# <u>28421</u>																																	
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1011</u> Instrument <u>1005</u> Mouth Alcohol Solution Lot # <u>2025-A</u> Acetone Stock Solution Lot # <u>2024-B</u>					
Simulator	Serial #	Lot #	Expiration																																	
0.000		N/A	N/A																																	
0.040																																				
0.100																																				
0.200																																				
0.300																																				
0.080 DGS	N/A																																			
☐ Post Calibration Adjustment Stability Checks			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6294</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>						Simulator	Serial Number	0.000	MP6294	Interferent	MP6290	0.050	MP5088	0.080	MP6292	0.200	MP6293																
Simulator	Serial Number																																			
0.000	MP6294																																			
Interferent	MP6290																																			
0.050	MP5088																																			
0.080	MP6292																																			
0.200	MP6293																																			
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☐ Other _____													
Simulator	Serial #	Lot #	Expiration																																	
0.050																																				
0.080																																				
0.200																																				
0.080 DGS	N/A																																			
Notes/Suggested Service: _____ _____ _____ _____ _____ _____			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																	
			<table border="0"><tr><td>Taylor Gutschow Tech Review / Date</td><td>Digitally signed by Taylor Gutschow Date: 2025.09.07 22:59:32 +0400</td><td>Shayla Platt Admin Review / Date</td><td>Digitally signed by Shayla Platt Date: 2025.09.11 08:24:48 -0400</td></tr></table>						Taylor Gutschow Tech Review / Date	Digitally signed by Taylor Gutschow Date: 2025.09.07 22:59:32 +0400	Shayla Platt Admin Review / Date	Digitally signed by Shayla Platt Date: 2025.09.11 08:24:48 -0400																								
Taylor Gutschow Tech Review / Date	Digitally signed by Taylor Gutschow Date: 2025.09.07 22:59:32 +0400	Shayla Platt Admin Review / Date	Digitally signed by Shayla Platt Date: 2025.09.11 08:24:48 -0400																																	

Stability Checks

80-005454

sat 9/3/25

0.050 g/210L	0.080 g/210L	0.200 g/210L	DGS 0.080 g/210L																																																																																																																																																
☒ 0.047 to 0.053 g/210L	☒ 0.077 to 0.083 g/210L	☒ 0.194 to 0.206 g/210L	☒ 0.077 to 0.083 g/210L <0.003 g/210L of Wet																																																																																																																																																
Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis																																																																																																																																																
OCALA PD Intoxilyzer - Alcohol Analyzer Model 8000 09/03/2025 Software: 8100.27 SN 80-005454 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:25</td></tr> <tr><td>Control Test</td><td>0.049</td><td>14:25</td></tr> <tr><td>Air Blank</td><td>0.003</td><td>14:26</td></tr> <tr><td>Control Test</td><td>0.050</td><td>14:27</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:27</td></tr> <tr><td>Control Test</td><td>0.043</td><td>14:28</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:28</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>2.0408</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:25	Control Test	0.049	14:25	Air Blank	0.003	14:26	Control Test	0.050	14:27	Air Blank	0.000	14:27	Control Test	0.043	14:28	Air Blank	0.000	14:28	Control Test Stats			Average	0.0490		Std Dev	0.0010		Rel Std Dev(%)	2.0408		OCALA PD Intoxilyzer - Alcohol Analyzer Model 8000 09/03/2025 Software: 8100.27 SN 80-005454 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:47</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:48</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:49</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:50</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:51</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:51</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0787</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:47	Control Test	0.078	14:48	Air Blank	0.000	14:49	Control Test	0.079	14:49	Air Blank	0.000	14:50	Control Test	0.079	14:51	Air Blank	0.000	14:51	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339		OCALA PD Intoxilyzer - Alcohol Analyzer Model 8000 09/03/2025 Software: 8100.27 SN 80-005454 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:42</td></tr> <tr><td>Control Test</td><td>0.198</td><td>14:43</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:43</td></tr> <tr><td>Control Test</td><td>0.198</td><td>14:44</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:44</td></tr> <tr><td>Control Test</td><td>0.199</td><td>14:45</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:45</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1983</td><td></td></tr> <tr><td>Std Dev</td><td>0.0005</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2911</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:42	Control Test	0.198	14:43	Air Blank	0.000	14:43	Control Test	0.198	14:44	Air Blank	0.000	14:44	Control Test	0.199	14:45	Air Blank	0.000	14:45	Control Test Stats			Average	0.1983		Std Dev	0.0005		Rel Std Dev(%)	0.2911		OCALA PD Intoxilyzer - Alcohol Analyzer Model 8000 09/03/2025 Software: 8100.27 SN 80-005454 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:18</td></tr> <tr><td>Control Test</td><td>0.077</td><td>14:18</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:19</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:20</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:21</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0777</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7434</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	14:18	Control Test	0.077	14:18	Air Blank	0.000	14:19	Control Test	0.078	14:19	Air Blank	0.000	14:20	Control Test	0.078	14:20	Air Blank	0.000	14:21	Control Test Stats			Average	0.0777		Std Dev	0.0006		Rel Std Dev(%)	0.7434	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:25																																																																																																																																																	
Control Test	0.049	14:25																																																																																																																																																	
Air Blank	0.003	14:26																																																																																																																																																	
Control Test	0.050	14:27																																																																																																																																																	
Air Blank	0.000	14:27																																																																																																																																																	
Control Test	0.043	14:28																																																																																																																																																	
Air Blank	0.000	14:28																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0490																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	2.0408																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:47																																																																																																																																																	
Control Test	0.078	14:48																																																																																																																																																	
Air Blank	0.000	14:49																																																																																																																																																	
Control Test	0.079	14:49																																																																																																																																																	
Air Blank	0.000	14:50																																																																																																																																																	
Control Test	0.079	14:51																																																																																																																																																	
Air Blank	0.000	14:51																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0787																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7339																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:42																																																																																																																																																	
Control Test	0.198	14:43																																																																																																																																																	
Air Blank	0.000	14:43																																																																																																																																																	
Control Test	0.198	14:44																																																																																																																																																	
Air Blank	0.000	14:44																																																																																																																																																	
Control Test	0.199	14:45																																																																																																																																																	
Air Blank	0.000	14:45																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1983																																																																																																																																																		
Std Dev	0.0005																																																																																																																																																		
Rel Std Dev(%)	0.2911																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:18																																																																																																																																																	
Control Test	0.077	14:18																																																																																																																																																	
Air Blank	0.000	14:19																																																																																																																																																	
Control Test	0.078	14:19																																																																																																																																																	
Air Blank	0.000	14:20																																																																																																																																																	
Control Test	0.078	14:20																																																																																																																																																	
Air Blank	0.000	14:21																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0777																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7434																																																																																																																																																		
 Operator's Signature	 Operator's Signature	 Operator's Signature	 Operator's Signature																																																																																																																																																

DGS

80-005454

DGS #2

Inadvertently ran DGS
stability twice — Both
tests were acceptable

SLT
9/3/2025

Ocala PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-005454
09/03/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	14:34
Control Test	0.078	14:34
Air Blank	0.000	14:34
Control Test	0.078	14:35
Air Blank	0.000	14:35
Control Test	0.078	14:36
Air Blank	0.000	14:36
Control Test Stats		
Average	0.0780	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: OCALA PD
Time of Inspection: 13:50

Date of Inspection: 09/04/2025

Serial Number: 80-005454
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.049	0.079	0.199	0.079
0.000	0.050	0.079	0.198	0.079
0.000	0.049	0.079	0.198	0.079
0.000	0.050	0.078	0.198	0.079
0.000	0.050	0.079	0.199	0.079
0.000	0.050	0.079	0.198	0.079
0.000	0.049	0.079	0.198	0.079
0.000	0.050	0.079	0.199	0.079
0.000	0.050	0.080	0.198	0.079
0.000	0.050	0.079	0.199	0.079

Standard Deviations	0.0004	0.0004	0.0005	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5



Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


LEANDRA HIGGINBOTHAM
Signature and Printed Name

09/04/2025
Date

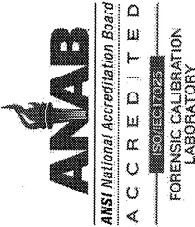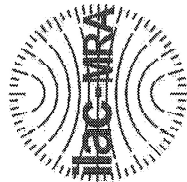

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-005454, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-005454</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>Ocala PD</u>	0.050 g/210 L	0.004
Calibration Date:	<u>09/04/2025</u>	0.080 g/210 L	0.004
Calibration Time:	<u>13:50</u>	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$)

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

09/04/2025

Date


LEANDRA HOGGINSBOTHAM,
Department Inspector

FDLE/ATP Form 69 October 2024
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality