

Alcohol Testing Program - Instrument Processing Sheet

Agency: PANAMA CITY BEACH PD Instrument Serial Number: 80-000838
 Date In: 1/6/2026 DI Completion Date: _____ Ship P/U H/D CMI EE

Intake By: <u>KTS</u> Date: <u>1/6/2026</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Dropped Off <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE <input type="checkbox"/> Training Instrument Visual Inspection <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/Accessories <input type="checkbox"/> Power Cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: Agency note. KTS 1/6/26	Quality Checks By: _____ Date: _____ <input type="checkbox"/> Breath Tube Screen <input type="checkbox"/> Replace External O-Rings <input type="checkbox"/> Instrument Set Up Verified <input type="checkbox"/> R-Value: _____ <input type="checkbox"/> Flow Verification (L/s) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547) <input type="checkbox"/> Barometric Pressure Check Gauge ID #: _____ Gauge: _____ Instrument: _____ <input type="checkbox"/> Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot#/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> </tr> </tbody> </table>	Simulator	Serial #	Lot#/Exp	0.050			0.080			0.200			0.080 DGS	N/A		Flow Adjustment By: _____ Flow Column #: _____ <input type="checkbox"/> 5L/min – 17mm <input type="checkbox"/> 15L/min – 53mm <input type="checkbox"/> 30L/min – 103mm <input type="checkbox"/> R-Value: _____ <input type="checkbox"/> Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)																					
Simulator	Serial #	Lot#/Exp																																				
0.050																																						
0.080																																						
0.200																																						
0.080 DGS	N/A																																					
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot#/Exp</th> <th>Maintenance</th> <th>By:</th> <th>Date:</th> </tr> </thead> <tbody> <tr> <td></td> <td></td> <td></td> <td><input type="checkbox"/> Battery Replacement</td> <td></td> <td></td> </tr> <tr> <td></td> <td></td> <td></td> <td><input type="checkbox"/> Dry Gas Regulator Replacement</td> <td></td> <td></td> </tr> <tr> <td></td> <td></td> <td></td> <td><input type="checkbox"/> Tank Sensor Tare</td> <td></td> <td></td> </tr> <tr> <td></td> <td></td> <td></td> <td><input type="checkbox"/> Breath Tube Replacement</td> <td></td> <td></td> </tr> <tr> <td></td> <td></td> <td></td> <td><input type="checkbox"/> Other:</td> <td></td> <td></td> </tr> </tbody> </table>			Simulator	Serial #	Lot#/Exp	Maintenance	By:	Date:				<input type="checkbox"/> Battery Replacement						<input type="checkbox"/> Dry Gas Regulator Replacement						<input type="checkbox"/> Tank Sensor Tare						<input type="checkbox"/> Breath Tube Replacement						<input type="checkbox"/> Other:		
Simulator	Serial #	Lot#/Exp	Maintenance	By:	Date:																																	
			<input type="checkbox"/> Battery Replacement																																			
			<input type="checkbox"/> Dry Gas Regulator Replacement																																			
			<input type="checkbox"/> Tank Sensor Tare																																			
			<input type="checkbox"/> Breath Tube Replacement																																			
			<input type="checkbox"/> Other:																																			
Optical Bench Adjustment By: _____		Department Inspection By: _____																																				
Barometric Pressure Gauge: _____ ID#: _____		Barometric Pressure ID#: _____																																				
<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge: _____ Instrument: _____ Mouth Alcohol Solution Lot #: _____ Exp: _____ Acetone Stock Solution Lot #: _____ Exp: _____								
Simulator	Serial #	Lot #	Expiration																																			
0.000		N/A	N/A																																			
0.040																																						
0.100																																						
0.200																																						
0.300																																						
0.080 DGS	N/A																																					
<input type="checkbox"/> Post Optical Bench Adjustment Stability Checks		<table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A																	
Simulator	Serial #	Lot #	Expiration																																			
0.050																																						
0.080																																						
0.200																																						
0.080 DGS	N/A																																					
Gauge ID #: _____ Gauge: _____ Instrument: _____		Attachments <input type="checkbox"/> Form 41 <input type="checkbox"/> Post-Stability Checks <input type="checkbox"/> Stability Checks <input type="checkbox"/> Flow Adjustment <input type="checkbox"/> Calibration Certificate <input type="checkbox"/> Form 40 <input type="checkbox"/> Optical Bench Adjustment <input type="checkbox"/> Other:																																				
Notes/Suggested Service: Corrected instrument settings. Previous inspection conducted on 12/3/25. KTS 1/6/26		<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC																																				
		<input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use																																				
		Shayla Platt <small>Digitally signed by Shayla Platt Date: 2026.01.07 13:20:30 -05'00'</small>	Shayla Platt <small>Digitally signed by Shayla Platt Date: 2026.01.07 13:20:43 -05'00'</small>																																			
		Tech Review	Admin Review																																			

80-000838

We have just received the instrument back from being repaired by CMI, the FDLE ATP. I was able to complete a valid monthly inspection on the instrument.

When attempting to use it for a breath test, the instrument goes directly into a Diagnostic Test, then requests a breath sample, without being able to input all the required arrest, observation, suspect, and arresting officer information.

A hard boot was conducted with the same results.

An upload was completed, and again, the same results.

If you could please complete the annual on the instrument since you have it, it would be greatly appreciated.

Bob Moss

Agency Inspector

Panama City Beach Police Department

850-625-0835

Bob.moss@pcbfl.gov

Instrument Processing Sheet

Agency: PANAMA CITY BEACH PD Instrument Serial Number: 80-000838
 Date In: 12/1/2025 DI Completion Date: 12/3/2025 Ship P/U H/D CMI EE

Intake By: <u>KTS</u> Date: <u>12/1/25</u>	Quality Checks By: <u>KTS</u> Date: <u>12/3/25</u>	Flow Adjustment By: _____ Date: _____															
<input checked="" type="checkbox"/> Annual <input type="checkbox"/> Registration <input checked="" type="checkbox"/> Return from CMI / EE <input type="checkbox"/> Return unworked <input type="checkbox"/> Training Visual Inspection <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/Accessories <input type="checkbox"/> Power Cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes:	<input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value: <u>117</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column #: <u>ATP103</u> 32 mm <u>0.160</u> (.139-.169) 36 mm <u>0.175</u> (.156-.190) 53 mm <u>0.250</u> (.228-.278) 103 mm <u>0.496</u> (.447-.547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1014</u> Instrument: <u>1015</u> <input checked="" type="checkbox"/> Stability Checks	Flow Column #: _____ <input type="checkbox"/> 5L/min – 17mm <input type="checkbox"/> 15L/min – 53mm <input type="checkbox"/> 30L/min – 103mm <input type="checkbox"/> R-Value: _____ <input type="checkbox"/> Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)															
	<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot#/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP5088</td> <td>202406K 6/19/2026</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> <td>202406L 6/19/2026</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>202406N 6/20/2026</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG510701 4/17/2027</td> </tr> </tbody> </table>	Simulator	Serial #	Lot#/Exp	0.050	MP5088	202406K 6/19/2026	0.080	MP5089	202406L 6/19/2026	0.200	MP5090	202406N 6/20/2026	0.080 DGS	N/A	AG510701 4/17/2027	Maintenance By: _____ Date: _____ <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement and Tank Sensor Tare <input type="checkbox"/> Breath Tube Replacement <input type="checkbox"/> Other:
Simulator	Serial #	Lot#/Exp															
0.050	MP5088	202406K 6/19/2026															
0.080	MP5089	202406L 6/19/2026															
0.200	MP5090	202406N 6/20/2026															
0.080 DGS	N/A	AG510701 4/17/2027															

Optical Bench Adjustment By: _____	Department Inspection By: <u>KTS</u>																																								
Barometric Pressure Gauge: _____ ID#: _____	Barometric Pressure ID#: <u>28427</u>																																								
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> <input type="checkbox"/> Post Optical Bench Adjustment Stability Checks	Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge: <u>1015</u> Instrument: <u>1016</u> Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>9/25/2027</u> Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>9/22/2027</u> <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP6289</td> </tr> <tr> <td>Interferent</td> <td>MP5087</td> </tr> <tr> <td>0.050</td> <td>MP5088</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> </tr> </tbody> </table>	Simulator	Serial Number	0.000	MP6289	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial #	Lot #	Expiration																																						
0.000		N/A	N/A																																						
0.040																																									
0.100																																									
0.200																																									
0.300																																									
0.080 DGS	N/A																																								
Simulator	Serial Number																																								
0.000	MP6289																																								
Interferent	MP5087																																								
0.050	MP5088																																								
0.080	MP5089																																								
0.200	MP5090																																								
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>	Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments <input checked="" type="checkbox"/> Form 41 <input type="checkbox"/> Post-Stability Checks <input checked="" type="checkbox"/> Stability Checks <input type="checkbox"/> Flow Adjustment <input checked="" type="checkbox"/> Calibration Certificate <input type="checkbox"/> Form 40 <input type="checkbox"/> Optical Bench Adjustment <input type="checkbox"/> Other:																				
Simulator	Serial #	Lot #	Expiration																																						
0.050																																									
0.080																																									
0.200																																									
0.080 DGS	N/A																																								
Barometric Pressure Gauge: _____ ID#: _____																																									

Notes/Suggested Service: Tech Review: Added checkmark to instrument return method. KTS 12/12/25	<input checked="" type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <hr/> <input checked="" type="checkbox"/> Return to/Place into Evidentiary Use <input type="checkbox"/> Remain Out of Evidentiary Use <hr/> <input checked="" type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use
	Digitally signed by Shayla Platt Shayla Platt Date: 2025.12.14 13:54:35 -05'00' <hr/> Digitally signed by Wen-Chi Pierson Wen-Chi Pierson Date: 2025.12.18 11:01:11 -05'00'
Tech Review	Admin Review

Stability Checks: Root Cause Analysis performed. Tightened simulator before re-conducting testing. KTS 12/3/25

PANAMA CITY BEACH PD
Intoxilyzer - Alcotest Analyzer
Model: 8000 SN 80-000838
12/03/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:18
Control Test	0.080	07:19
Air Blank	0.000	07:19
Control Test	0.080	07:19
Air Blank	0.000	07:20
Control Test	0.081	07:20
Air Blank	0.000	07:21
Control Test Status		
Average	0.0803	
Std Dev	0.0065	
Rel Std Dev(%)	0.7137	

Katrina S
Operator's Signature

Dry


PANAMA CITY BEACH PD
Intoxilyzer - Alcotest Analyzer
Model: 8000 SN 80-000838
12/03/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:24
Control Test	0.040	07:25
Air Blank	0.000	07:25
Control Test	0.049	07:26
Air Blank	0.000	07:26
Control Test	0.048	07:27
Air Blank	0.000	07:27
Control Test Status		
Average	0.0457	
Std Dev	0.0049	
Rel Std Dev(%)	10.8019	

Katrina S
Operator's Signature

PANAMA CITY BEACH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000838
12/03/2025
Software: 8100.27


Test	g/210L	Time
Air Blank	0.000	07:29
Control Test	0.048	07:29
Air Blank	0.000	07:30
Control Test	0.048	07:31
Air Blank	0.000	07:31
Control Test	0.049	07:32
Air Blank	0.000	07:32
Control Test Stats		
Average	0.0483	
Std Dev	0.0006	
Rel Std Dev(%)	1.1945	

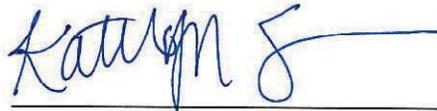

Operator's Signature

PANAMA CITY BEACH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000838
12/03/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:34
Control Test	0.078	07:34
Air Blank	0.000	07:35
Control Test	0.078	07:36
Air Blank	0.000	07:36
Control Test	0.079	07:37
Air Blank	0.000	07:37
Control Test Stats		
Average	0.0783	
Std Dev	0.0006	
Rel Std Dev(%)	0.7370	

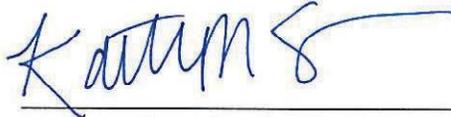


Operator's Signature

Wet

PANAMA CITY BEACH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000838
12/03/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:38
Control Test	0.199	07:39
Air Blank	0.000	07:40
Control Test	0.199	07:40
Air Blank	0.000	07:41
Control Test	0.199	07:42
Air Blank	0.000	07:42
Control Test Stats		
Average	0.1990	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

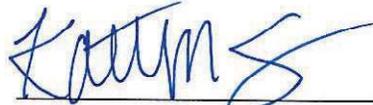

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PANAMA CITY BEACH PD
Time of Inspection: 10:08

Date of Inspection: 12/03/2025

Serial Number: 80-000838
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG510701 Exp: 04/17/2027
0.000	0.048	0.078	0.199	0.080
0.000	0.048	0.078	0.198	0.080
0.000	0.048	0.078	0.199	0.081
0.000	0.048	0.078	0.198	0.080
0.000	0.048	0.078	0.199	0.081
0.000	0.048	0.078	0.198	0.080
0.000	0.048	0.078	0.199	0.081
0.000	0.048	0.078	0.198	0.080
0.000	0.048	0.078	0.199	0.081
0.000	0.048	0.078	0.198	0.080
0.000	0.048	0.078	0.199	0.081

Standard Deviations	0.0000	0.0000	0.0005	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0002 Number of Simulators Used: 5

Remarks:

The above instrument complies () does not comply () with Chapter 11D-8, FAC.



I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.



KATIE T SPEARIN

Signature and Printed Name

12/03/2025
Date

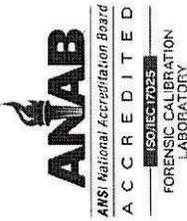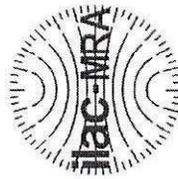

Florida Department of Law Enforcement
 Alcohol Testing Program
 2331 Phillips Road
 Tallahassee, FL 32308

Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000838, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000838</u>	UNCERTAINTY* ±	
Owning Agency:	<u>PANAMA CITY BEACH PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>12/03/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>10:08</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.
 Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
 *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).
 The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
 This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Katie Spearin

12/03/2025

Date
KATIE TSPEARIN,
 Department Inspector

FDLE/ATP Form 69 October 2024
 Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Stability Checks

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L ≤0.003 g/210L of Wet																																																																																																																																																																		
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>																																																																																																																																																																		
Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis																																																																																																																																																																		
<p>PANAMA CITY BEACH PD Intoxilyzer - Alcohol Analyzer Model 8100 SN 60-001838 07/16/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:36</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:36</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:37</td></tr> <tr><td>Control Test</td><td>0.048</td><td>13:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:38</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:39</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:39</td></tr> <tr><td>Control Test</td><td>0.000</td><td>13:39</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0497</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel. Std Dev.(%)</td><td>1.1663</td><td></td></tr> </tbody> </table> <p><i>Katz</i> Operator's Signature</p>	Test	g/210L	Time	Air Blank	0.000	13:36	Control Test	0.049	13:36	Air Blank	0.000	13:37	Control Test	0.048	13:37	Air Blank	0.000	13:38	Control Test	0.049	13:39	Air Blank	0.000	13:39	Control Test	0.000	13:39	Control Test Stats			Average	0.0497		Std Dev	0.0006		Rel. Std Dev.(%)	1.1663		<p>PANAMA CITY BEACH PD Intoxilyzer - Alcohol Analyzer Model 8100 SN 60-001838 07/16/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:41</td></tr> <tr><td>Control Test</td><td>0.077</td><td>13:42</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:43</td></tr> <tr><td>Control Test</td><td>0.077</td><td>13:43</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:44</td></tr> <tr><td>Control Test</td><td>0.077</td><td>13:44</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:45</td></tr> <tr><td>Control Test</td><td>0.000</td><td>13:45</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0770</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel. Std Dev.(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> <p><i>Katz</i> Operator's Signature</p>	Test	g/210L	Time	Air Blank	0.000	13:41	Control Test	0.077	13:42	Air Blank	0.000	13:43	Control Test	0.077	13:43	Air Blank	0.000	13:44	Control Test	0.077	13:44	Air Blank	0.000	13:45	Control Test	0.000	13:45	Control Test Stats			Average	0.0770		Std Dev	0.0000		Rel. Std Dev.(%)	0.0000		<p>PANAMA CITY BEACH PD Intoxilyzer - Alcohol Analyzer Model 8100 SN 60-001838 07/16/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:40</td></tr> <tr><td>Control Test</td><td>0.196</td><td>13:40</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:47</td></tr> <tr><td>Control Test</td><td>0.196</td><td>13:47</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:48</td></tr> <tr><td>Control Test</td><td>0.195</td><td>13:48</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:49</td></tr> <tr><td>Control Test</td><td>0.195</td><td>13:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:50</td></tr> <tr><td>Control Test</td><td>0.000</td><td>13:50</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.1957</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel. Std Dev.(%)</td><td>0.2891</td><td></td></tr> </tbody> </table> <p><i>Katz</i> Operator's Signature</p>	Test	g/210L	Time	Air Blank	0.000	13:40	Control Test	0.196	13:40	Air Blank	0.000	13:47	Control Test	0.196	13:47	Air Blank	0.000	13:48	Control Test	0.195	13:48	Air Blank	0.000	13:49	Control Test	0.195	13:49	Air Blank	0.000	13:50	Control Test	0.000	13:50	Control Test Stats			Average	0.1957		Std Dev	0.0006		Rel. Std Dev.(%)	0.2891		<p>PANAMA CITY BEACH PD Intoxilyzer - Alcohol Analyzer Model 8100 SN 60-001838 07/16/2025 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:29</td></tr> <tr><td>Control Test</td><td>0.080</td><td>13:30</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:30</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:30</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:31</td></tr> <tr><td>Control Test</td><td>0.080</td><td>13:31</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:32</td></tr> <tr><td>Control Test</td><td>0.000</td><td>13:32</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0797</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel. Std Dev.(%)</td><td>0.7247</td><td></td></tr> </tbody> </table> <p><i>Katz</i> Operator's Signature</p>	Test	g/210L	Time	Air Blank	0.000	13:29	Control Test	0.080	13:30	Air Blank	0.000	13:30	Control Test	0.079	13:30	Air Blank	0.000	13:31	Control Test	0.080	13:31	Air Blank	0.000	13:32	Control Test	0.000	13:32	Control Test Stats			Average	0.0797		Std Dev	0.0006		Rel. Std Dev.(%)	0.7247	
Test	g/210L	Time																																																																																																																																																																			
Air Blank	0.000	13:36																																																																																																																																																																			
Control Test	0.049	13:36																																																																																																																																																																			
Air Blank	0.000	13:37																																																																																																																																																																			
Control Test	0.048	13:37																																																																																																																																																																			
Air Blank	0.000	13:38																																																																																																																																																																			
Control Test	0.049	13:39																																																																																																																																																																			
Air Blank	0.000	13:39																																																																																																																																																																			
Control Test	0.000	13:39																																																																																																																																																																			
Control Test Stats																																																																																																																																																																					
Average	0.0497																																																																																																																																																																				
Std Dev	0.0006																																																																																																																																																																				
Rel. Std Dev.(%)	1.1663																																																																																																																																																																				
Test	g/210L	Time																																																																																																																																																																			
Air Blank	0.000	13:41																																																																																																																																																																			
Control Test	0.077	13:42																																																																																																																																																																			
Air Blank	0.000	13:43																																																																																																																																																																			
Control Test	0.077	13:43																																																																																																																																																																			
Air Blank	0.000	13:44																																																																																																																																																																			
Control Test	0.077	13:44																																																																																																																																																																			
Air Blank	0.000	13:45																																																																																																																																																																			
Control Test	0.000	13:45																																																																																																																																																																			
Control Test Stats																																																																																																																																																																					
Average	0.0770																																																																																																																																																																				
Std Dev	0.0000																																																																																																																																																																				
Rel. Std Dev.(%)	0.0000																																																																																																																																																																				
Test	g/210L	Time																																																																																																																																																																			
Air Blank	0.000	13:40																																																																																																																																																																			
Control Test	0.196	13:40																																																																																																																																																																			
Air Blank	0.000	13:47																																																																																																																																																																			
Control Test	0.196	13:47																																																																																																																																																																			
Air Blank	0.000	13:48																																																																																																																																																																			
Control Test	0.195	13:48																																																																																																																																																																			
Air Blank	0.000	13:49																																																																																																																																																																			
Control Test	0.195	13:49																																																																																																																																																																			
Air Blank	0.000	13:50																																																																																																																																																																			
Control Test	0.000	13:50																																																																																																																																																																			
Control Test Stats																																																																																																																																																																					
Average	0.1957																																																																																																																																																																				
Std Dev	0.0006																																																																																																																																																																				
Rel. Std Dev.(%)	0.2891																																																																																																																																																																				
Test	g/210L	Time																																																																																																																																																																			
Air Blank	0.000	13:29																																																																																																																																																																			
Control Test	0.080	13:30																																																																																																																																																																			
Air Blank	0.000	13:30																																																																																																																																																																			
Control Test	0.079	13:30																																																																																																																																																																			
Air Blank	0.000	13:31																																																																																																																																																																			
Control Test	0.080	13:31																																																																																																																																																																			
Air Blank	0.000	13:32																																																																																																																																																																			
Control Test	0.000	13:32																																																																																																																																																																			
Control Test Stats																																																																																																																																																																					
Average	0.0797																																																																																																																																																																				
Std Dev	0.0006																																																																																																																																																																				
Rel. Std Dev.(%)	0.7247																																																																																																																																																																				

PANAMA CITY BEACH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000838
07/16/2025
Software: 8100.27

Flow Rate Calibration*****

1: Rate (Liters/min) = 5

SQRT(Diff)) = 6.242

2: Rate (Liters/min) = 15

SQRT(Diff)) = 11.703

3: Rate (Liters/min) = 30

SQRT(Diff)) = 23.598

Dependent Data Scale Factor = 100000 L/min

Independent Data Scale Factor = 256

Rounded Slope = 551

Rounded Intercept = -287813

Correlation = 0.99541

#1

SUA

7/16/25

PANAMA CITY BEACH PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000838
07/16/2025
Software: 8100.27

Flow Rate Calibration*****

- 1: Rate (Liters/min) = 5
SQRT(Diff)) = 6.164
- 2: Rate (Liters/min) = 15
SQRT(Diff)) = 11.312
- 3: Rate (Liters/min) = 30
SQRT(Diff)) = 21.539

Dependent Data Scale Factor = 100000 L/min

Independent Data Scale Factor = 256

Rounded Slope = 626

Rounded Intercept = -418733

Correlation = 0.99730

#2

SUA

7/16/25

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: PANAMA CITY BEACH PD
Time of Inspection: 14:36

Date of Inspection: 07/17/2025

Serial Number: 80-000838
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.048	0.078	0.195	0.079
0.000	0.048	0.078	0.195	0.079
0.000	0.048	0.079	0.195	0.079
0.000	0.049	0.079	0.195	0.079
0.000	0.048	0.079	0.195	0.078
0.000	0.049	0.079	0.195	0.079
0.000	0.048	0.078	0.195	0.078
0.000	0.049	0.079	0.195	0.079
0.000	0.048	0.078	0.195	0.079
0.000	0.048	0.079	0.195	0.079


Standard Deviations	0.0004	0.0005	0.0000	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies () does not comply () with Chapter 11D-8, FAC.



I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

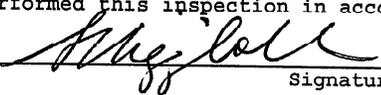

LEANDRA HIGGINBOTHAM

Signature and Printed Name

07/17/2025
Date

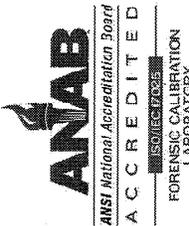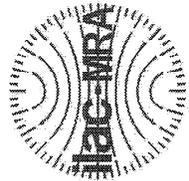

Florida Department of Law Enforcement
 Alcohol Testing Program
 2331 Phillips Road
 Tallahassee, FL 32308

Calibration Certificate

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000838, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-000838</u>	UNCERTAINTY* ±	
Owning Agency:	<u>PANAMA CITY BEACH PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>07/17/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>14:36</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
 *Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.
 This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.


 LEANDRA HIGGINBOTHAM,
 Department Inspector

07/17/2025
 Date

FDLE/ATP Form 69 October 2024
 Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality