

INSTRUMENT PROCESSING SHEET

Agency Miami Dade SO

S/N 80-006895

Florida Department of
Law Enforcement

Date In 9/8/2025

DI Completion Date 09/16/2025

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By <u>SLH</u> Date <u>9/8/2025</u> ☐ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ 	Quality Checks By <u>SLH</u> Date <u>9/15/2025</u> ☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>220</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP103</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.234</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28421</u> ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> <tr> <td>0.050</td> <td>MP6291</td> <td>202406K 06/19/2026</td> </tr> <tr> <td>0.080</td> <td>MP6292</td> <td>202406L 06/19/2026</td> </tr> <tr> <td>0.200</td> <td>MP6293</td> <td>202406N 06/20/2026</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG429602 10/22/2026</td> </tr> </table>	Simulator	Serial #	Lot #/Exp	0.050	MP6291	202406K 06/19/2026	0.080	MP6292	202406L 06/19/2026	0.200	MP6293	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	Flow Calibration By _____ Date _____ Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By <u>SLH</u> Date <u>9/9/2025</u> ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☒ Other <u>CMI loaded forms remotely.</u> <u>Form 40 was generated to verify.</u> _____ _____ _____ 																																														
Simulator	Serial #	Lot #/Exp																																																													
0.050	MP6291	202406K 06/19/2026																																																													
0.080	MP6292	202406L 06/19/2026																																																													
0.200	MP6293	202406N 06/20/2026																																																													
0.080 DGS	N/A	AG429602 10/22/2026																																																													
Calibration Adjustment By _____ Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table> Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ 		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection By <u>SLH</u> Barometric Pressure ID# <u>28421</u> Gauge <u>1015</u> Instrument <u>1015</u> Mouth Alcohol Solution Lot # <u>2025-A</u> Acetone Stock Solution Lot # <u>2024-B</u> <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> <tr> <td>0.000</td> <td>MP6289</td> </tr> <tr> <td>Interferent</td> <td>MP5087</td> </tr> <tr> <td>0.050</td> <td>MP6291</td> </tr> <tr> <td>0.080</td> <td>MP6292</td> </tr> <tr> <td>0.200</td> <td>MP6293</td> </tr> </table> Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☒ Form 40 ☒ Other <u>Note to load form</u> ☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.09.18 16:25:10 -0400 Phil Nicodemo Digitally signed by Phil Nicodemo Date: 2025.09.23 15:06:07 -0400 Tech Review / Date _____ Admin Review / Date _____		Simulator	Serial Number	0.000	MP6289	Interferent	MP5087	0.050	MP6291	0.080	MP6292	0.200	MP6293
Simulator	Serial #	Lot #	Expiration																																																												
0.000		N/A	N/A																																																												
0.040																																																															
0.100																																																															
0.200																																																															
0.300																																																															
0.080 DGS	N/A																																																														
Simulator	Serial #	Lot #	Expiration																																																												
0.050																																																															
0.080																																																															
0.200																																																															
0.080 DGS	N/A																																																														
Simulator	Serial Number																																																														
0.000	MP6289																																																														
Interferent	MP5087																																																														
0.050	MP6291																																																														
0.080	MP6292																																																														
0.200	MP6293																																																														

Taylor

80-006895 need 5 forms
installed.

You've already done the
MDSO switch on it

- Rick.

Rec'd 9/8/25
SMT

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI DADE SO
Time of Inspection: 12:34

Date of Inspection: 09/09/2025

Serial Number: 80-006895
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08g/210L Test (g/210L) Lot#:_____ Exp:_____	0.20g/210L Test (g/210L) Lot#:_____ Exp:_____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:_____ Exp:_____

Number of Simulators Used: _____

Remarks:

BYPASS AI TO VERIFY FORMS LOADED SUCCESSFULLY

✓ *SUB*

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

LEANDRA HIGGINBOTHAM

LEANDRA HIGGINBOTHAM

Signature and Printed Name

09/09/2025
Date

80-006895 sub-
9.15.25

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L ≤0.003 g/210L of Wet																																																																																																																																																
✓	✓	✓	✓																																																																																																																																																
Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis																																																																																																																																																
MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 09/15/2025 Software: 8100.27 SN 80-006895 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>11:43</td></tr><tr><td>Control Test</td><td>0.048</td><td>11:44</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:44</td></tr><tr><td>Control Test</td><td>0.049</td><td>11:45</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:46</td></tr><tr><td>Control Test</td><td>0.049</td><td>11:46</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:47</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0487</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.1863</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	11:43	Control Test	0.048	11:44	Air Blank	0.000	11:44	Control Test	0.049	11:45	Air Blank	0.000	11:46	Control Test	0.049	11:46	Air Blank	0.000	11:47	Control Test Stats			Average	0.0487		Std Dev	0.0006		Rel Std Dev(%)	1.1863		MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 09/15/2025 Software: 8100.27 SN 80-006895 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>11:54</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:55</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:55</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:56</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:56</td></tr><tr><td>Control Test</td><td>0.080</td><td>11:57</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:58</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0793</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.7277</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	11:54	Control Test	0.079	11:55	Air Blank	0.000	11:55	Control Test	0.079	11:56	Air Blank	0.000	11:56	Control Test	0.080	11:57	Air Blank	0.000	11:58	Control Test Stats			Average	0.0793		Std Dev	0.0006		Rel Std Dev(%)	0.7277		MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 09/15/2025 Software: 8100.27 SN 80-006895 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>11:59</td></tr><tr><td>Control Test</td><td>0.198</td><td>12:00</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:01</td></tr><tr><td>Control Test</td><td>0.199</td><td>12:01</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:02</td></tr><tr><td>Control Test</td><td>0.198</td><td>12:03</td></tr><tr><td>Air Blank</td><td>0.000</td><td>12:03</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.1983</td><td></td></tr><tr><td>Std Dev</td><td>0.0006</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>0.2911</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	11:59	Control Test	0.198	12:00	Air Blank	0.000	12:01	Control Test	0.199	12:01	Air Blank	0.000	12:02	Control Test	0.198	12:03	Air Blank	0.000	12:03	Control Test Stats			Average	0.1983		Std Dev	0.0006		Rel Std Dev(%)	0.2911		MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 09/15/2025 Software: 8100.27 SN 80-006895 <table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.000</td><td>11:48</td></tr><tr><td>Control Test</td><td>0.081</td><td>11:49</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:49</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:50</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:50</td></tr><tr><td>Control Test</td><td>0.079</td><td>11:50</td></tr><tr><td>Air Blank</td><td>0.000</td><td>11:51</td></tr><tr><td>Control Test Stats</td><td></td><td></td></tr><tr><td>Average</td><td>0.0797</td><td></td></tr><tr><td>Std Dev</td><td>0.0012</td><td></td></tr><tr><td>Rel Std Dev(%)</td><td>1.4494</td><td></td></tr></tbody></table>	Test	g/210L	Time	Air Blank	0.000	11:48	Control Test	0.081	11:49	Air Blank	0.000	11:49	Control Test	0.079	11:50	Air Blank	0.000	11:50	Control Test	0.079	11:50	Air Blank	0.000	11:51	Control Test Stats			Average	0.0797		Std Dev	0.0012		Rel Std Dev(%)	1.4494	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:43																																																																																																																																																	
Control Test	0.048	11:44																																																																																																																																																	
Air Blank	0.000	11:44																																																																																																																																																	
Control Test	0.049	11:45																																																																																																																																																	
Air Blank	0.000	11:46																																																																																																																																																	
Control Test	0.049	11:46																																																																																																																																																	
Air Blank	0.000	11:47																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0487																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	1.1863																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:54																																																																																																																																																	
Control Test	0.079	11:55																																																																																																																																																	
Air Blank	0.000	11:55																																																																																																																																																	
Control Test	0.079	11:56																																																																																																																																																	
Air Blank	0.000	11:56																																																																																																																																																	
Control Test	0.080	11:57																																																																																																																																																	
Air Blank	0.000	11:58																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0793																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7277																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:59																																																																																																																																																	
Control Test	0.198	12:00																																																																																																																																																	
Air Blank	0.000	12:01																																																																																																																																																	
Control Test	0.199	12:01																																																																																																																																																	
Air Blank	0.000	12:02																																																																																																																																																	
Control Test	0.198	12:03																																																																																																																																																	
Air Blank	0.000	12:03																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1983																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2911																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:48																																																																																																																																																	
Control Test	0.081	11:49																																																																																																																																																	
Air Blank	0.000	11:49																																																																																																																																																	
Control Test	0.079	11:50																																																																																																																																																	
Air Blank	0.000	11:50																																																																																																																																																	
Control Test	0.079	11:50																																																																																																																																																	
Air Blank	0.000	11:51																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0797																																																																																																																																																		
Std Dev	0.0012																																																																																																																																																		
Rel Std Dev(%)	1.4494																																																																																																																																																		
Miguel Operator's Signature	Miguel Operator's Signature	Miguel Operator's Signature	Miguel Operator's Signature																																																																																																																																																

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI DADE SO
Time of Inspection: 18:29

Date of Inspection: 09/16/2025

Serial Number: 80-006895
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: 202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#: 202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#: 202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: AG429602 Exp: 10/22/2026
0.000	0.049	0.078	0.196	0.081
0.000	0.048	0.078	0.196	0.082
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.049	0.080	0.196	0.081
0.000	0.049	0.079	0.197	0.081
0.000	0.048	0.079	0.197	0.081


Standard Deviations	0.0004	0.0006	0.0004	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

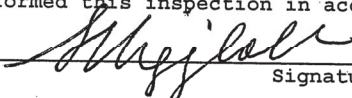
LEANDRA HIGGINBOTHAM
Signature and Printed Name

09/16/2025
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006895, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-006895	UNCERTAINTY* $\pm$	
Owning Agency:	MIAMI DADE SO	0.050 g/210 L	0.004
Calibration Date:	09/16/2025	0.080 g/210 L	0.004
Calibration Time:	18:29	0.200 g/210 L	0.007
		0.080 g/210 L Dry Gas Control	0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

09/16/2025
Date

LEANDRA HIGGINBOTHAM,
Department Inspector

FDLE/ATP Form 69 October 2024
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

INSTRUMENT PROCESSING SHEET

Agency Miami Dade Sheriff's Office

S/N 80-006895

Florida Department of
Law Enforcement

Date In 03/14/2025 DI Completion Date 04/23/2025

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Date 4/22/2025	Quality Checks	By TDG	Date 4/22/2025																																																												
☒ Annual ☒ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value 215 ☒ Flow Verification (L/s) Flow Column # ATP101 32 mm 0.156 (.139 - .169) 36 mm 0.171 (.156 - .190) 53 mm 0.242 (.228 - .278) 103 mm 0.503 (.447 - .547) ☒ Barometric Pressure Check Gauge ID # 28199 ☒ Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP6286</td> <td>202406K 06/19/2026</td> </tr> <tr> <td>0.080</td> <td>MP6287</td> <td>202406L 06/19/2026</td> </tr> <tr> <td>0.200</td> <td>MP6288</td> <td>202406N 06/20/2026</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG429602 10/22/2026</td> </tr> </tbody> </table>			Simulator	Serial #	Lot #/Exp	0.050	MP6286	202406K 06/19/2026	0.080	MP6287	202406L 06/19/2026	0.200	MP6288	202406N 06/20/2026	0.080 DGS	N/A	AG429602 10/22/2026	Flow Calibration By _____ Date _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547) Maintenance By _____ Date _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Breath Tube Replacement ☐ Other _____ _____ _____ _____																																												
Simulator	Serial #	Lot #/Exp																																																															
0.050	MP6286	202406K 06/19/2026																																																															
0.080	MP6287	202406L 06/19/2026																																																															
0.200	MP6288	202406N 06/20/2026																																																															
0.080 DGS	N/A	AG429602 10/22/2026																																																															
Calibration Adjustment Barometric Pressure Gauge _____ ID # _____ <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table> ☐ Post Calibration Adjustment Stability Checks <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Department Inspection Barometric Pressure ID# 28199 Gauge 1020 Instrument 1020 Mouth Alcohol Solution Lot # 2024-A Acetone Stock Solution Lot # 2023-B <table border="1" style="width:100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP6284</td> </tr> <tr> <td>Interferent</td> <td>MP6285</td> </tr> <tr> <td>0.050</td> <td>MP6286</td> </tr> <tr> <td>0.080</td> <td>MP6287</td> </tr> <tr> <td>0.200</td> <td>MP6288</td> </tr> </tbody> </table>			Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288
Simulator	Serial #	Lot #	Expiration																																																														
0.000		N/A	N/A																																																														
0.040																																																																	
0.100																																																																	
0.200																																																																	
0.300																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial #	Lot #	Expiration																																																														
0.050																																																																	
0.080																																																																	
0.200																																																																	
0.080 DGS	N/A																																																																
Simulator	Serial Number																																																																
0.000	MP6284																																																																
Interferent	MP6285																																																																
0.050	MP6286																																																																
0.080	MP6287																																																																
0.200	MP6288																																																																
Notes/Suggested Service: After the 2025 inspection, FDLE was notified the agency's name had changed to Miami Dade Sheriff's Office. Retained the instrument to update the registration. (TDG 4/22/25) _____ _____ _____ _____			Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☒ Other Form 51																																																														
Notes/Suggested Service: After the 2025 inspection, FDLE was notified the agency's name had changed to Miami Dade Sheriff's Office. Retained the instrument to update the registration. (TDG 4/22/25)			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Digitally signed by Phil Nicodemo Date: 2025.05.02 12:54:58 -0400 Phil Nicodemo Shayla Platt Digitally signed by Shayla Platt Date: 2025.05.03 10:04:00 Admin Review Date																																																														
Tech Review / Date _____			Admin Review / Date _____																																																														



Florida Department of
Law Enforcement

REQUEST FOR REGISTRATION

MAKE AND MODEL OF INSTRUMENT: Intoxilyzer 8000

SERIAL NUMBER: 80-006895

OWNING AGENCY: Miami Dade Sheriff's Office

DATE OF DEPARTMENT INSPECTION: 04/23/2025

AGENCY INSPECTOR: Richard Closius

ADDRESS: 1519 NW 79th Avenue

CITY, STATE, ZIP: Miami, FL 33126

TELEPHONE NUMBER: 305-458-3991

FAX NUMBER: N/A

EMAIL ADDRESS (if available): rlclosius@mdso.com

For Program Office Use Only:

- ☒ Registration Issued
- ☒ Instrument Added to Evidentiary Instrument Database
- ☐ Instrument Added to Monthly Statistics Database
- ☒ Contact Information Added to Instrument Database

Stability Checks

0.05g/210L 0.047 to 0.053 ✓	0.08g/210L 0.077 to 0.083 ✓	0.20g/210L 0.194 to 0.206 ✓	DGS 0.08g/210L 0.077 to 0.083 ✓ ≤0.003 of Wet ✓																																																																																																																																																
MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006895 04/22/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:04</td></tr> <tr><td>Control Test</td><td>0.049</td><td>15:05</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:05</td></tr> <tr><td>Control Test</td><td>0.048</td><td>15:06</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:06</td></tr> <tr><td>Control Test</td><td>0.048</td><td>15:07</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:08</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0483</td><td></td></tr> <tr><td>Std Dev</td><td>0.0016</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1945</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	15:04	Control Test	0.049	15:05	Air Blank	0.000	15:05	Control Test	0.048	15:06	Air Blank	0.000	15:06	Control Test	0.048	15:07	Air Blank	0.000	15:08	Control Test Stats			Average	0.0483		Std Dev	0.0016		Rel Std Dev(%)	1.1945		MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006895 04/22/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:11</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:12</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:12</td></tr> <tr><td>Control Test</td><td>0.078</td><td>15:13</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:13</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:14</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:15</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0787</td><td></td></tr> <tr><td>Std Dev</td><td>0.0016</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	15:11	Control Test	0.079	15:12	Air Blank	0.000	15:12	Control Test	0.078	15:13	Air Blank	0.000	15:13	Control Test	0.079	15:14	Air Blank	0.000	15:15	Control Test Stats			Average	0.0787		Std Dev	0.0016		Rel Std Dev(%)	0.7339		MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006895 04/22/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:55</td></tr> <tr><td>Control Test</td><td>0.197</td><td>14:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:57</td></tr> <tr><td>Control Test</td><td>0.197</td><td>14:57</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:58</td></tr> <tr><td>Control Test</td><td>0.196</td><td>14:59</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:59</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.1967</td><td></td></tr> <tr><td>Std Dev</td><td>0.0016</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2936</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	14:55	Control Test	0.197	14:56	Air Blank	0.000	14:57	Control Test	0.197	14:57	Air Blank	0.000	14:58	Control Test	0.196	14:59	Air Blank	0.000	14:59	Control Test Stats			Average	0.1967		Std Dev	0.0016		Rel Std Dev(%)	0.2936		MIAMI DADE SO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006895 04/22/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:16</td></tr> <tr><td>Control Test</td><td>0.081</td><td>15:16</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:16</td></tr> <tr><td>Control Test</td><td>0.080</td><td>15:17</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:17</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:18</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:18</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0800</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.2500</td><td></td></tr> </tbody> </table> Operator's Signature 	Test	g/210L	Time	Air Blank	0.000	15:16	Control Test	0.081	15:16	Air Blank	0.000	15:16	Control Test	0.080	15:17	Air Blank	0.000	15:17	Control Test	0.079	15:18	Air Blank	0.000	15:18	Control Test Stats			Average	0.0800		Std Dev	0.0010		Rel Std Dev(%)	1.2500	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:04																																																																																																																																																	
Control Test	0.049	15:05																																																																																																																																																	
Air Blank	0.000	15:05																																																																																																																																																	
Control Test	0.048	15:06																																																																																																																																																	
Air Blank	0.000	15:06																																																																																																																																																	
Control Test	0.048	15:07																																																																																																																																																	
Air Blank	0.000	15:08																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0483																																																																																																																																																		
Std Dev	0.0016																																																																																																																																																		
Rel Std Dev(%)	1.1945																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:11																																																																																																																																																	
Control Test	0.079	15:12																																																																																																																																																	
Air Blank	0.000	15:12																																																																																																																																																	
Control Test	0.078	15:13																																																																																																																																																	
Air Blank	0.000	15:13																																																																																																																																																	
Control Test	0.079	15:14																																																																																																																																																	
Air Blank	0.000	15:15																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0787																																																																																																																																																		
Std Dev	0.0016																																																																																																																																																		
Rel Std Dev(%)	0.7339																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:55																																																																																																																																																	
Control Test	0.197	14:56																																																																																																																																																	
Air Blank	0.000	14:57																																																																																																																																																	
Control Test	0.197	14:57																																																																																																																																																	
Air Blank	0.000	14:58																																																																																																																																																	
Control Test	0.196	14:59																																																																																																																																																	
Air Blank	0.000	14:59																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1967																																																																																																																																																		
Std Dev	0.0016																																																																																																																																																		
Rel Std Dev(%)	0.2936																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:16																																																																																																																																																	
Control Test	0.081	15:16																																																																																																																																																	
Air Blank	0.000	15:16																																																																																																																																																	
Control Test	0.080	15:17																																																																																																																																																	
Air Blank	0.000	15:17																																																																																																																																																	
Control Test	0.079	15:18																																																																																																																																																	
Air Blank	0.000	15:18																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0800																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	1.2500																																																																																																																																																		

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI DADE SO
Time of Inspection: 12:11

Date of Inspection: 04/23/2025

Serial Number: 80-006895
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: 202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#: 202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#: 202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: AG429602 Exp: 10/22/2026
0.000	0.048	0.078	0.196	0.082
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.048	0.078	0.197	0.081
0.000	0.048	0.079	0.196	0.081
0.000	0.048	0.079	0.197	0.081
0.000	0.049	0.079	0.197	0.081


Standard Deviations	0.0003	0.0005	0.0004	0.0003
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

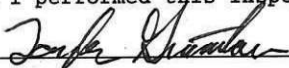

Signature and Printed Name

TAYLOR D GUTSCHOW

04/23/2025
Date

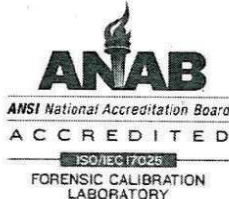

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006895, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006895</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>MIAMI DADE SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>04/23/2025</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>12:11</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

04/23/2025

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality



INSTRUMENT PROCESSING SHEET

Agency Miami Dade PDS/N 80-006895Florida Department of
Law EnforcementDate In 03/14/2025 DI Completion Date 03/25/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By TDG	Date	Quality Checks	By TDG	Date	Flow Calibration	By	Date																												
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____		03/25/2025	☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>215</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP101</u> 32 mm <u>0.160</u> (.139 - .169) 36 mm <u>0.175</u> (.156 - .190) 53 mm <u>0.250</u> (.228 - .278) 103 mm <u>0.511</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>26932</u> ☒ Stability Checks		03/25/2025	Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																														
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP6286</td><td>202303K 03/29/2025</td></tr><tr><td>0.080</td><td>MP6287</td><td>202303L 03/29/2025</td></tr><tr><td>0.200</td><td>MP6288</td><td>202304C 04/05/2025</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG429602 10/22/2026</td></tr></tbody></table>	Simulator	Serial #	Lot #/Exp	0.050	MP6286	202303K 03/29/2025	0.080	MP6287	202303L 03/29/2025	0.200	MP6288	202304C 04/05/2025	0.080 DGS	N/A	AG429602 10/22/2026																		
Simulator	Serial #	Lot #/Exp																																		
0.050	MP6286	202303K 03/29/2025																																		
0.080	MP6287	202303L 03/29/2025																																		
0.200	MP6288	202304C 04/05/2025																																		
0.080 DGS	N/A	AG429602 10/22/2026																																		
Calibration Adjustment			Department Inspection																																	
By _____			By <u>TDG</u>																																	
Barometric Pressure Gauge ID # _____			Barometric Pressure ID# <u>26932</u>																																	
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1018</u> Instrument <u>1018</u> Mouth Alcohol Solution Lot # <u>2024-A</u> Acetone Stock Solution Lot # <u>2023-B</u>					
Simulator	Serial #	Lot #	Expiration																																	
0.000		N/A	N/A																																	
0.040																																				
0.100																																				
0.200																																				
0.300																																				
0.080 DGS	N/A																																			
☐ Post Calibration Adjustment Stability Checks			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6284</td></tr><tr><td>Interferent</td><td>MP6285</td></tr><tr><td>0.050</td><td>MP6286</td></tr><tr><td>0.080</td><td>MP6287</td></tr><tr><td>0.200</td><td>MP6288</td></tr></tbody></table>						Simulator	Serial Number	0.000	MP6284	Interferent	MP6285	0.050	MP6286	0.080	MP6287	0.200	MP6288																
Simulator	Serial Number																																			
0.000	MP6284																																			
Interferent	MP6285																																			
0.050	MP6286																																			
0.080	MP6287																																			
0.200	MP6288																																			
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Calibration Adjustment ☐ Post-Stability Checks ☐ Flow Calibration ☐ Form 40 ☐ Other _____													
Simulator	Serial #	Lot #	Expiration																																	
0.050																																				
0.080																																				
0.200																																				
0.080 DGS	N/A																																			
Notes/Suggested Service: _____ _____ _____ _____ _____ _____ _____ _____ _____ _____			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Destinee Armstrong Digitally signed by Destinee Armstrong Date: 2025.03.26 14:48:51 -0400 Shayla Platt Digitally signed by Shayla Platt Date: 2025.03.27 08:43:33 -0400 Tech Review / Date _____ Admin Review / Date _____																																	

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																																
MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model 8000 03/25/2025 Software: 8100.27 SN 80-006895 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:11</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:12</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:13</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:13</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:14</td></tr> <tr><td>Control Test</td><td>0.049</td><td>13:15</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:15</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> Operator's Signature <i>M</i>	Test	g/210L	Time	Air Blank	0.000	13:11	Control Test	0.049	13:12	Air Blank	0.000	13:13	Control Test	0.049	13:13	Air Blank	0.000	13:14	Control Test	0.049	13:15	Air Blank	0.000	13:15	Control Test Stats			Average	0.0490		Std Dev	0.0000		Rel Std Dev(%)	0.0000		MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model 8000 03/25/2025 Software: 8100.27 SN 80-006895 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:19</td></tr> <tr><td>Control Test</td><td>0.080</td><td>13:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:20</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:21</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:22</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:22</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:23</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0793</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7277</td><td></td></tr> </tbody> </table> Operator's Signature <i>M</i>	Test	g/210L	Time	Air Blank	0.000	13:19	Control Test	0.080	13:20	Air Blank	0.000	13:20	Control Test	0.079	13:21	Air Blank	0.000	13:22	Control Test	0.079	13:22	Air Blank	0.000	13:23	Control Test Stats			Average	0.0793		Std Dev	0.0006		Rel Std Dev(%)	0.7277		MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model 8000 03/25/2025 Software: 8100.27 SN 80-006895 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:58</td></tr> <tr><td>Control Test</td><td>0.196</td><td>12:59</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:00</td></tr> <tr><td>Control Test</td><td>0.196</td><td>13:00</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:01</td></tr> <tr><td>Control Test</td><td>0.197</td><td>13:02</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:02</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1963</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2941</td><td></td></tr> </tbody> </table> Operator's Signature <i>M</i>	Test	g/210L	Time	Air Blank	0.000	12:58	Control Test	0.196	12:59	Air Blank	0.000	13:00	Control Test	0.196	13:00	Air Blank	0.000	13:01	Control Test	0.197	13:02	Air Blank	0.000	13:02	Control Test Stats			Average	0.1963		Std Dev	0.0006		Rel Std Dev(%)	0.2941		MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model 8000 03/25/2025 Software: 8100.27 SN 80-006895 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>12:11</td></tr> <tr><td>Control Test</td><td>0.081</td><td>12:11</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:12</td></tr> <tr><td>Control Test</td><td>0.080</td><td>12:12</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:13</td></tr> <tr><td>Control Test</td><td>0.080</td><td>12:13</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:14</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0803</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7187</td><td></td></tr> </tbody> </table> Operator's Signature <i>M</i>	Test	g/210L	Time	Air Blank	0.000	12:11	Control Test	0.081	12:11	Air Blank	0.000	12:12	Control Test	0.080	12:12	Air Blank	0.000	12:13	Control Test	0.080	12:13	Air Blank	0.000	12:14	Control Test Stats			Average	0.0803		Std Dev	0.0006		Rel Std Dev(%)	0.7187	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	13:11																																																																																																																																																	
Control Test	0.049	13:12																																																																																																																																																	
Air Blank	0.000	13:13																																																																																																																																																	
Control Test	0.049	13:13																																																																																																																																																	
Air Blank	0.000	13:14																																																																																																																																																	
Control Test	0.049	13:15																																																																																																																																																	
Air Blank	0.000	13:15																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0490																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	13:19																																																																																																																																																	
Control Test	0.080	13:20																																																																																																																																																	
Air Blank	0.000	13:20																																																																																																																																																	
Control Test	0.079	13:21																																																																																																																																																	
Air Blank	0.000	13:22																																																																																																																																																	
Control Test	0.079	13:22																																																																																																																																																	
Air Blank	0.000	13:23																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0793																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7277																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:58																																																																																																																																																	
Control Test	0.196	12:59																																																																																																																																																	
Air Blank	0.000	13:00																																																																																																																																																	
Control Test	0.196	13:00																																																																																																																																																	
Air Blank	0.000	13:01																																																																																																																																																	
Control Test	0.197	13:02																																																																																																																																																	
Air Blank	0.000	13:02																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1963																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2941																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	12:11																																																																																																																																																	
Control Test	0.081	12:11																																																																																																																																																	
Air Blank	0.000	12:12																																																																																																																																																	
Control Test	0.080	12:12																																																																																																																																																	
Air Blank	0.000	12:13																																																																																																																																																	
Control Test	0.080	12:13																																																																																																																																																	
Air Blank	0.000	12:14																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0803																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7187																																																																																																																																																		

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: MIAMI DADE PD
Time of Inspection: 15:07

Date of Inspection: 03/25/2025

Serial Number: 80-006895
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202303K Exp: 03/29/2025	0.08g/210L Test (g/210L) Lot#:202303L Exp: 03/29/2025	0.20g/210L Test (g/210L) Lot#:202304C Exp: 04/05/2025	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG429602 Exp: 10/22/2026
0.000	0.050	0.079	0.196	0.082
0.000	0.049	0.079	0.196	0.081
0.000	0.050	0.079	0.196	0.081
0.000	0.050	0.079	0.197	0.082
0.000	0.050	0.079	0.196	0.081
0.000	0.049	0.079	0.196	0.082
0.000	0.050	0.079	0.196	0.081
0.000	0.049	0.079	0.196	0.081
0.000	0.051	0.079	0.196	0.081
0.000	0.050	0.079	0.196	0.081


Standard Deviations	0.0006	0.0000	0.0003	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.


I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

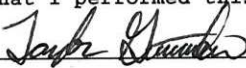

TAYLOR D GUTSCHOW

Signature and Printed Name

03/25/2025
Date

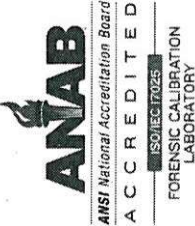

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
4700 Terminal Drive, Suite 1
Ft. Myers, FL 33907

This is to certify the calibration of Intoxilyzer 8000 serial number 80-006895, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-006895</u>	UNCERTAINTY* $\pm$
Owning Agency:	<u>MIAMI DADE PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>03/25/2025</u>	0.080 g/ 210 L 0.004
Calibration Time:	<u>15:07</u>	0.200 g/ 210 L 0.007
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

03/25/2025

Date

TAYLOR D GUTSCHOW,
Department Inspector

FDLE/ATP Form 69 December 2021

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality