

Alcohol Testing Program - Instrument Processing Sheet

Agency: LAKE ALFRED PD

Instrument Serial Number: 80-001048

Date In: 5/20/2026

DI Completion Date: 5/29/26

☒ Ship

☐ P/U

☐ H/D

☐ CMI

☐ EE

Intake By: <u>WKP</u> Date: <u>5/20/2026</u>		Quality Checks By: <u>IAS</u> Date: <u>5/22/2026</u>		Flow Adjustment By: _____																																
☒ Annual ☐ Dropped Off ☐ Registration ☒ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☐ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>280</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP105</u> 32 mm <u>.156</u> (.139-.169) 36 mm <u>.171</u> (.156-.190) 53 mm <u>.246</u> (.228-.278) 103 mm <u>.519</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1011</u> Instrument: <u>1016</u> ☒ Stability Checks		Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)																																
		<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 20%;">Serial #</th> <th style="width: 25%;">Lot#/Exp</th> </tr> </thead> <tbody> <tr> <td rowspan="2">0.050</td> <td rowspan="2">MP5088</td> <td>202603A</td> </tr> <tr> <td>3/3/2028</td> </tr> <tr> <td rowspan="2">0.080</td> <td rowspan="2">MP5089</td> <td>202603B</td> </tr> <tr> <td>3/3/2028</td> </tr> <tr> <td rowspan="2">0.200</td> <td rowspan="2">MP5090</td> <td>202603D</td> </tr> <tr> <td>3/4/2028</td> </tr> <tr> <td rowspan="2">0.080 DGS</td> <td rowspan="2">N/A</td> <td>23725080A1</td> </tr> <tr> <td>10/5/2027</td> </tr> </tbody> </table>		Simulator	Serial #	Lot#/Exp	0.050	MP5088	202603A	3/3/2028	0.080	MP5089	202603B	3/3/2028	0.200	MP5090	202603D	3/4/2028	0.080 DGS	N/A	23725080A1	10/5/2027	<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 60%;">Maintenance</th> <th style="width: 40%;">By: _____ Date: _____</th> </tr> </thead> <tbody> <tr> <td>☐ Battery Replacement</td> <td></td> </tr> <tr> <td>☐ Dry Gas Regulator Replacement</td> <td></td> </tr> <tr> <td>☐ Tank Sensor Tare</td> <td></td> </tr> <tr> <td>☐ Breath Tube Replacement</td> <td></td> </tr> <tr> <td>☐ Other:</td> <td></td> </tr> </tbody> </table>		Maintenance	By: _____ Date: _____	☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Tank Sensor Tare		☐ Breath Tube Replacement		☐ Other:	
Simulator	Serial #	Lot#/Exp																																		
0.050	MP5088	202603A																																		
		3/3/2028																																		
0.080	MP5089	202603B																																		
		3/3/2028																																		
0.200	MP5090	202603D																																		
		3/4/2028																																		
0.080 DGS	N/A	23725080A1																																		
		10/5/2027																																		
Maintenance	By: _____ Date: _____																																			
☐ Battery Replacement																																				
☐ Dry Gas Regulator Replacement																																				
☐ Tank Sensor Tare																																				
☐ Breath Tube Replacement																																				
☐ Other:																																				
Optical Bench Adjustment By: _____		Department Inspection By: <u>KTS</u>																																		
Barometric Pressure Gauge: _____ ID#: _____		Barometric Pressure ID#: <u>28427</u>																																		
Simulator	Serial #	Lot #	Expiration	Gauge: <u>1010</u> Instrument: <u>1014</u>																																
0.000		N/A	N/A	Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>9/25/27</u>																																
0.040				Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>9/22/27</u>																																
0.100				Simulator	Serial Number																															
0.200				0.000	MP5086																															
0.300				Interferent	MP5087																															
0.080 DGS	N/A			0.050	MP5088																															
☐ Post Optical Bench Adjustment Stability Checks				0.080	MP5089																															
Simulator	Serial #	Lot #	Expiration	0.200	MP5090																															
0.050				Attachments																																
0.080				☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Optical Bench Adjustment																																
0.200																																				
0.080 DGS	N/A																																			
Gauge ID #: _____				☐ Post-Stability Checks ☐ Flow Adjustment ☒ Form 40 ☐ Other:																																
Gauge: _____ Instrument: _____																																				
Notes/Suggested Service:				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																
				<table style="width: 100%;"> <tr> <td style="width: 50%;"> LeAndra Higginbotham Digitally signed by LeAndra Higginbotham Date: 2026.06.04 15:52:40 -04'00' </td> <td style="width: 50%;"> Wen-Chi Pierson Digitally signed by Wen-Chi Pierson Date: 2026.06.05 09:02:48 -04'00' </td> </tr> </table>		LeAndra Higginbotham Digitally signed by LeAndra Higginbotham Date: 2026.06.04 15:52:40 -04'00'	Wen-Chi Pierson Digitally signed by Wen-Chi Pierson Date: 2026.06.05 09:02:48 -04'00'																													
LeAndra Higginbotham Digitally signed by LeAndra Higginbotham Date: 2026.06.04 15:52:40 -04'00'	Wen-Chi Pierson Digitally signed by Wen-Chi Pierson Date: 2026.06.05 09:02:48 -04'00'																																			
				<table style="width: 100%;"> <tr> <td style="width: 50%; text-align: center;">Tech Review</td> <td style="width: 50%; text-align: center;">Admin Review</td> </tr> </table>		Tech Review	Admin Review																													
Tech Review	Admin Review																																			

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: LAKE ALFRED PD
Time of Inspection: 09:54

Date of Inspection: 05/22/2026

Serial Number: 80-001048
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:

Number of Simulators Used: _____

Remarks:

BYPASS AI FOR OPERATION, COMPLIANCE UNDETERMINED.

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Isabella Smykowski

ISABELLA A SMYKOWSKI

Signature and Printed Name

05/22/2026
Date

Stability Checks

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L ≤0.003 g/210L of Wet																																																																																																																																																						
LAKE ALFRED PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001048 05/22/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:34</td></tr> <tr><td>Control Test</td><td>0.051</td><td>13:35</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:35</td></tr> <tr><td>Control Test</td><td>0.050</td><td>13:36</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:37</td></tr> <tr><td>Control Test</td><td>0.051</td><td>13:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:38</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0507</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>1.1395</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	13:34	Control Test	0.051	13:35	Air Blank	0.000	13:35	Control Test	0.050	13:36	Air Blank	0.000	13:37	Control Test	0.051	13:37	Air Blank	0.000	13:38	Control Test Stats			Average	0.0507		Std Dev	0.0006		Rel Std Dev(%)	1.1395		LAKE ALFRED PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001048 05/22/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:23</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:24</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:25</td></tr> <tr><td>Control Test</td><td>0.078</td><td>13:25</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:26</td></tr> <tr><td>Control Test</td><td>0.079</td><td>13:27</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:27</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0787</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	13:23	Control Test	0.079	13:24	Air Blank	0.000	13:25	Control Test	0.078	13:25	Air Blank	0.000	13:26	Control Test	0.079	13:27	Air Blank	0.000	13:27	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339		LAKE ALFRED PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001048 05/22/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>13:18</td></tr> <tr><td>Control Test</td><td>0.198</td><td>13:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:19</td></tr> <tr><td>Control Test</td><td>0.198</td><td>13:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:20</td></tr> <tr><td>Control Test</td><td>0.198</td><td>13:21</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:21</td></tr> <tr><td>Control Test</td><td>0.198</td><td>13:21</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>13:22</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.1980</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	13:18	Control Test	0.198	13:19	Air Blank	0.000	13:19	Control Test	0.198	13:19	Air Blank	0.000	13:20	Control Test	0.198	13:21	Air Blank	0.000	13:21	Control Test	0.198	13:21	Air Blank	0.000	13:22	Control Test Stats			Average	0.1980		Std Dev	0.0000		Rel Std Dev(%)	0.0000		LAKE ALFRED PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-001048 05/22/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:05</td></tr> <tr><td>Control Test</td><td>0.079</td><td>14:06</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:06</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:06</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:07</td></tr> <tr><td>Control Test</td><td>0.078</td><td>14:07</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:08</td></tr> <tr><td colspan="3">Control Test Stats</td></tr> <tr><td>Average</td><td>0.0783</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7370</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	14:05	Control Test	0.079	14:06	Air Blank	0.000	14:06	Control Test	0.078	14:06	Air Blank	0.000	14:07	Control Test	0.078	14:07	Air Blank	0.000	14:08	Control Test Stats			Average	0.0783		Std Dev	0.0006		Rel Std Dev(%)	0.7370	
Test	g/210L	Time																																																																																																																																																							
Air Blank	0.000	13:34																																																																																																																																																							
Control Test	0.051	13:35																																																																																																																																																							
Air Blank	0.000	13:35																																																																																																																																																							
Control Test	0.050	13:36																																																																																																																																																							
Air Blank	0.000	13:37																																																																																																																																																							
Control Test	0.051	13:37																																																																																																																																																							
Air Blank	0.000	13:38																																																																																																																																																							
Control Test Stats																																																																																																																																																									
Average	0.0507																																																																																																																																																								
Std Dev	0.0006																																																																																																																																																								
Rel Std Dev(%)	1.1395																																																																																																																																																								
Test	g/210L	Time																																																																																																																																																							
Air Blank	0.000	13:23																																																																																																																																																							
Control Test	0.079	13:24																																																																																																																																																							
Air Blank	0.000	13:25																																																																																																																																																							
Control Test	0.078	13:25																																																																																																																																																							
Air Blank	0.000	13:26																																																																																																																																																							
Control Test	0.079	13:27																																																																																																																																																							
Air Blank	0.000	13:27																																																																																																																																																							
Control Test Stats																																																																																																																																																									
Average	0.0787																																																																																																																																																								
Std Dev	0.0006																																																																																																																																																								
Rel Std Dev(%)	0.7339																																																																																																																																																								
Test	g/210L	Time																																																																																																																																																							
Air Blank	0.000	13:18																																																																																																																																																							
Control Test	0.198	13:19																																																																																																																																																							
Air Blank	0.000	13:19																																																																																																																																																							
Control Test	0.198	13:19																																																																																																																																																							
Air Blank	0.000	13:20																																																																																																																																																							
Control Test	0.198	13:21																																																																																																																																																							
Air Blank	0.000	13:21																																																																																																																																																							
Control Test	0.198	13:21																																																																																																																																																							
Air Blank	0.000	13:22																																																																																																																																																							
Control Test Stats																																																																																																																																																									
Average	0.1980																																																																																																																																																								
Std Dev	0.0000																																																																																																																																																								
Rel Std Dev(%)	0.0000																																																																																																																																																								
Test	g/210L	Time																																																																																																																																																							
Air Blank	0.000	14:05																																																																																																																																																							
Control Test	0.079	14:06																																																																																																																																																							
Air Blank	0.000	14:06																																																																																																																																																							
Control Test	0.078	14:06																																																																																																																																																							
Air Blank	0.000	14:07																																																																																																																																																							
Control Test	0.078	14:07																																																																																																																																																							
Air Blank	0.000	14:08																																																																																																																																																							
Control Test Stats																																																																																																																																																									
Average	0.0783																																																																																																																																																								
Std Dev	0.0006																																																																																																																																																								
Rel Std Dev(%)	0.7370																																																																																																																																																								

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: LAKE ALFRED PD

Time of Inspection: 13:21

Date of Inspection: 05/29/2026

Serial Number: 80-001048

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) <u>23825090A1</u> <i>RTS</i> Lot#: <u>23835080A1</u> Exp: 10/05/2027
0.000	0.050	0.079	0.198	0.080
0.000	0.049	0.079	0.198	0.080
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.080	0.199	0.080
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.079	0.198	0.080
0.000	0.050	0.079	0.198	0.080
0.000	0.050	0.079	0.198	0.080
0.000	0.050	0.080	0.198	0.080
0.000	0.050	0.079	0.199	0.080


Standard Deviations	0.0003	0.0004	0.0005	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

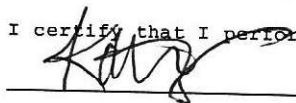

KATIE T SPEARIN
Signature and Printed Name

05/29/2026
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001048, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number: 80-001048

Owning Agency: LAKE ALFRED PD

Calibration Date: 05/29/2026

Calibration Time: 13:21

UNCERTAINTY* $\pm$

0.050 g/ 210 L 0.004

0.080 g/ 210 L 0.004

0.200 g/ 210 L 0.008

0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/29/2026

Date

KATIE T SPEARIN,
Department Inspector