

Alcohol Testing Program - Instrument Processing Sheet

Agency: NEW PORT RICHEY PD

Instrument Serial Number: 80-001064

Date In: 6/11/2026

DI Completion Date: 6/17/26

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By: <u>SLH</u> Date: <u>6/11/2026</u>		Quality Checks By: <u>KTS</u> Date: <u>6/12/2026</u>		Flow Adjustment By: _____																																	
☒ Annual ☐ Dropped Off ☐ Registration ☐ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☐ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>192</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP103</u> 32 mm <u>0.144</u> (.139-.169) 36 mm <u>0.160</u> (.156-.190) 53 mm <u>0.230</u> (.228-.278) 103 mm <u>0.492</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1012</u> Instrument: <u>1015</u> ☒ Stability Checks		Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)																																	
						<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot#/Exp</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td>MP5088</td> <td>202603A 3/3/2028</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> <td>202603B 3/3/2028</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>202603D 3/4/2028</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>23825080A1 10/5/2027</td> </tr> </tbody> </table>		Simulator	Serial #	Lot#/Exp	0.050	MP5088	202603A 3/3/2028	0.080	MP5089	202603B 3/3/2028	0.200	MP5090	202603D 3/4/2028	0.080 DGS	N/A	23825080A1 10/5/2027	Maintenance By: _____ Date: _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Tank Sensor Tare ☐ Breath Tube Replacement ☐ Other:														
						Simulator	Serial #	Lot#/Exp																													
						0.050	MP5088	202603A 3/3/2028																													
						0.080	MP5089	202603B 3/3/2028																													
0.200	MP5090	202603D 3/4/2028																																			
0.080 DGS	N/A	23825080A1 10/5/2027																																			
Optical Bench Adjustment By: _____		Department Inspection By: <u>KTS</u>																																			
Barometric Pressure Gauge: _____ ID#: _____		Barometric Pressure ID#: <u>28427</u>																																			
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge: <u>1009</u> Instrument: <u>1012</u> Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>9/25/27</u> Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>9/22/27</u>							
Simulator	Serial #	Lot #	Expiration																																		
0.000		N/A	N/A																																		
0.040																																					
0.100																																					
0.200																																					
0.300																																					
0.080 DGS	N/A																																				
☐ Post Optical Bench Adjustment Stability Checks		<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> </thead> <tbody> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </tbody> </table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> </thead> <tbody> <tr> <td>0.000</td> <td>MP6289</td> </tr> <tr> <td>Interferent</td> <td>MP5087</td> </tr> <tr> <td>0.050</td> <td>MP5088</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> </tr> </tbody> </table>		Simulator	Serial Number	0.000	MP6289	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial #	Lot #	Expiration																																		
0.050																																					
0.080																																					
0.200																																					
0.080 DGS	N/A																																				
Simulator	Serial Number																																				
0.000	MP6289																																				
Interferent	MP5087																																				
0.050	MP5088																																				
0.080	MP5089																																				
0.200	MP5090																																				
Gauge ID #: _____ Gauge: _____ Instrument: _____		Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Adjustment ☒ Calibration Certificate ☐ Form 40 ☐ Optical Bench Adjustment ☐ Other:																																			
Notes/Suggested Service:		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC																																			
		☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use																																			
		☒ Conduct an Agency Inspection Before Evidentiary Use																																			
Wen-Chi Pierson		LeAndra Higginbotham																																			
Tech Review		Admin Review																																			

Stability Checks

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L ≤0.003 g/210L of Wet
NEW PORT RICHEY PD Intoxilyzer - Alcohol Analyzer Model 8000 06/12/2026 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 08:45 Control Test 0.049 08:46 Air Blank 0.000 08:46 Control Test 0.049 08:47 Air Blank 0.000 08:47 Control Test 0.050 08:48 Air Blank 0.000 08:49 Control Test Stats Average 0.0493 Std Dev 0.0006 Rel Std Dev(%) 1.1703 Operator's Signature 	NEW PORT RICHEY PD Intoxilyzer - Alcohol Analyzer Model 8000 06/12/2026 Software: 8100.27 SN 60-001064 Test g/210L Time Air Blank 0.000 08:50 Control Test 0.078 08:51 Air Blank 0.000 08:51 Control Test 0.078 08:52 Air Blank 0.000 08:52 Control Test 0.078 08:53 Air Blank 0.000 08:54 Control Test Stats Average 0.0780 Std Dev 0.0000 Rel Std Dev(%) 0.0000 Operator's Signature 	NEW PORT RICHEY PD Intoxilyzer - Alcohol Analyzer Model 8000 06/12/2026 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 08:55 Control Test 0.197 08:56 Air Blank 0.000 08:57 Control Test 0.196 08:57 Air Blank 0.000 08:58 Control Test 0.196 08:59 Air Blank 0.000 08:59 Control Test Stats Average 0.1963 Std Dev 0.0006 Rel Std Dev(%) 0.2941 Operator's Signature 	NEW PORT RICHEY PD Intoxilyzer - Alcohol Analyzer Model 8000 06/12/2026 Software: 8100.27 SN 80-001064 Test g/210L Time Air Blank 0.000 08:38 Control Test 0.080 08:39 Air Blank 0.000 08:39 Control Test 0.081 08:40 Air Blank 0.000 08:40 Control Test 0.081 08:40 Air Blank 0.000 08:41 Control Test Stats Average 0.0807 Std Dev 0.0006 Rel Std Dev(%) 0.7157 Operator's Signature 

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: NEW PORT RICHEY PD
Time of Inspection: 12:38

Date of Inspection: 06/17/2026

Serial Number: 80-001064
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

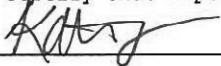
Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202603D Exp: 03/04/2028	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:23825080A1 Exp: 10/05/2027
0.000	0.050	0.078	0.196	0.081
0.000	0.050	0.078	0.197	0.081
0.000	0.049	0.078	0.197	0.080
0.000	0.050	0.078	0.197	0.080
0.000	0.050	0.078	0.197	0.080
0.000	0.050	0.078	0.198	0.081
0.000	0.049	0.078	0.198	0.080
0.000	0.050	0.078	0.197	0.081
0.000	0.050	0.078	0.198	0.081
0.000	0.050	0.078	0.198	0.080
Standard Deviations	0.0004	0.0000	0.0006	0.0005

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

 KATIE T SPEARIN
Signature and Printed Name

06/17/2026
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001064 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001064</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>NEW PORT RICHEY PD</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>06/17/2026</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>12:38</u>	0.200 g/ 210 L	0.008
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Date

06/17/2026

KATIE T SPEARIN,
Department Inspector

FDLE/ATP Form 69 January 2026
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality