

Alcohol Testing Program - Instrument Processing Sheet

Agency: Orange County SO

Instrument Serial Number: 80-001255

Date In: 4/3/2026

DI Completion Date: 05/18/2026

☒ Ship

☐ P/U

☐ H/D

☐ CMI

☐ EE

Intake By: <u>IAS</u> Date: <u>4/3/2026</u>		Quality Checks By: <u>SLH</u> Date: <u>5/8/2026</u>		Flow Adjustment By: _____																																							
☒ Annual ☐ Dropped Off ☐ Registration ☒ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☐ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>154</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP 105</u> 32 mm <u>0.152</u> (.139-.169) 36 mm <u>0.171</u> (.156-.190) 53 mm <u>0.250</u> (.228-.278) 103 mm <u>0.519</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1010</u> Instrument: <u>1013</u> ☒ Stability Checks		Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm (.139-.169) 36 mm (.156-.190) 53 mm (.228-.278) 103 mm (.447-.547)																																							
						<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 25%;">Serial #</th> <th style="width: 20%;">Lot#/Exp</th> <th style="width: 40%;">Maintenance By: <u>SLH</u> Date: <u>5/6/26</u></th> </tr> <tr> <td rowspan="2">0.050</td> <td rowspan="2">MP5088</td> <td>202603A</td> <td>☒ Battery Replacement</td> </tr> <tr> <td>03/03/2028</td> <td>☐ Dry Gas Regulator Replacement</td> </tr> <tr> <td rowspan="2">0.080</td> <td rowspan="2">MP5089</td> <td>202603B</td> <td>☐ Tank Sensor Tare</td> </tr> <tr> <td>03/03/2028</td> <td>☐ Breath Tube Replacement</td> </tr> <tr> <td rowspan="2">0.200</td> <td rowspan="2">MP5090</td> <td>202603D</td> <td>☒ Other:</td> </tr> <tr> <td>03/04/2028</td> <td>CMI remoted into instrument 4/7 and 4/8/26 to reset PW to default. Post battery change PW needed to be reset again. CMI was successful with the reset on 5/7/2026. SLH 5/8/26</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>23825080A1</td> <td></td> </tr> <tr> <td></td> <td></td> <td>10/05/2027</td> <td></td> </tr> </table>		Simulator	Serial #	Lot#/Exp	Maintenance By: <u>SLH</u> Date: <u>5/6/26</u>	0.050	MP5088	202603A	☒ Battery Replacement	03/03/2028	☐ Dry Gas Regulator Replacement	0.080	MP5089	202603B	☐ Tank Sensor Tare	03/03/2028	☐ Breath Tube Replacement	0.200	MP5090	202603D	☒ Other:	03/04/2028	CMI remoted into instrument 4/7 and 4/8/26 to reset PW to default. Post battery change PW needed to be reset again. CMI was successful with the reset on 5/7/2026. SLH 5/8/26	0.080 DGS	N/A	23825080A1				10/05/2027							
						Simulator	Serial #	Lot#/Exp	Maintenance By: <u>SLH</u> Date: <u>5/6/26</u>																																		
						0.050	MP5088	202603A	☒ Battery Replacement																																		
								03/03/2028	☐ Dry Gas Regulator Replacement																																		
0.080	MP5089	202603B	☐ Tank Sensor Tare																																								
		03/03/2028	☐ Breath Tube Replacement																																								
0.200	MP5090	202603D	☒ Other:																																								
		03/04/2028	CMI remoted into instrument 4/7 and 4/8/26 to reset PW to default. Post battery change PW needed to be reset again. CMI was successful with the reset on 5/7/2026. SLH 5/8/26																																								
0.080 DGS	N/A	23825080A1																																									
		10/05/2027																																									
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 25%;">Serial #</th> <th style="width: 20%;">Lot #</th> <th style="width: 40%;">Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <th style="width: 50%;">Simulator</th> <th style="width: 50%;">Serial Number</th> </tr> <tr> <td>0.000</td> <td>MP5086</td> </tr> <tr> <td>Interferent</td> <td>MP5087</td> </tr> <tr> <td>0.050</td> <td>MP5088</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> </tr> </table>		Simulator	Serial Number	0.000	MP5086	Interferent	MP5087	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial #	Lot #	Expiration																																								
0.000		N/A	N/A																																								
0.040																																											
0.100																																											
0.200																																											
0.300																																											
0.080 DGS	N/A																																										
Simulator	Serial Number																																										
0.000	MP5086																																										
Interferent	MP5087																																										
0.050	MP5088																																										
0.080	MP5089																																										
0.200	MP5090																																										
☐ Post Optical Bench Adjustment Stability Checks		Gauge: <u>1019</u> Instrument: <u>1021</u> Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>09/25/2027</u> Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>09/22/2027</u>																																									
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <th style="width: 15%;">Simulator</th> <th style="width: 25%;">Serial #</th> <th style="width: 20%;">Lot #</th> <th style="width: 40%;">Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <th colspan="2">Attachments</th> </tr> <tr> <td style="width: 50%; vertical-align: top;"> ☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Optical Bench Adjustment </td> <td style="width: 50%; vertical-align: top;"> ☐ Post-Stability Checks ☐ Flow Adjustment ☒ Form 40 ☐ Other: </td> </tr> </table>		Attachments		☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Optical Bench Adjustment	☐ Post-Stability Checks ☐ Flow Adjustment ☒ Form 40 ☐ Other:																
Simulator	Serial #	Lot #	Expiration																																								
0.050																																											
0.080																																											
0.200																																											
0.080 DGS	N/A																																										
Attachments																																											
☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Optical Bench Adjustment	☐ Post-Stability Checks ☐ Flow Adjustment ☒ Form 40 ☐ Other:																																										
Gauge ID #: _____ Gauge: _____ Instrument: _____																																											
Notes/Suggested Service: Tech Review: 1) Added date to 'Other' note in the Maintenance portion of the IPS. 2) Additional note regarding CMI remote service in the Maintenance portion of the IPS: 4/30/26 discovered needed forms loaded also and by 5/7/26 forms had been loaded, PW reset and cleared error. 3) On 6/1/26, spoke to Agency Inspector about PW and was able to demonstrate Level 2 and 3 PW operated as expected. PW=password; SLH 6/1/2026		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																									
		<table border="0" style="width: 100%;"> <tr> <td style="width: 50%; text-align: center;"> Digitally signed by Taylor Gutschow Date: 2026.06.03 15:27:20 -04'00' </td> <td style="width: 50%; text-align: center;"> Digitally signed by Wen-Chi Pierson Date: 2026.06.04 14:57:33 -04'00' </td> </tr> </table>		Digitally signed by Taylor Gutschow Date: 2026.06.03 15:27:20 -04'00'	Digitally signed by Wen-Chi Pierson Date: 2026.06.04 14:57:33 -04'00'																																						
		Digitally signed by Taylor Gutschow Date: 2026.06.03 15:27:20 -04'00'	Digitally signed by Wen-Chi Pierson Date: 2026.06.04 14:57:33 -04'00'																																								
<table border="0" style="width: 100%;"> <tr> <td style="width: 50%; text-align: center;"> Taylor Gutschow Tech Review </td> <td style="width: 50%; text-align: center;"> Wen-Chi Pierson Admin Review </td> </tr> </table>		Taylor Gutschow Tech Review	Wen-Chi Pierson Admin Review																																								
Taylor Gutschow Tech Review	Wen-Chi Pierson Admin Review																																										

Florida Department of Law Enforcement
Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: ✖

Serial Number: 80-001255

Time of Inspection: 12:23

Date of Inspection: 05/08/2026

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot #: Exp:

Number of Simulators Used: _____

Remarks:

AI BYPASS FOR OPERATION, COMPLIANCE UNDETERMINED

** Orange County SO - printed prior to verifying^{agency name} setup menu. Added after ward. sub*

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

Mijjala

LEANDRA HIGGINBOTHAM

Signature and Printed Name

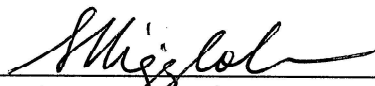
05/08/2026

Date

- noted

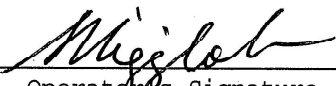
ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001255
05/08/2026
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:49
Control Test	0.048	12:50
Air Blank	0.000	12:50
Control Test	0.048	12:51
Air Blank	0.000	12:52
Control Test	0.049	12:52
Air Blank	0.000	12:53
Control Test Stats		
Average	0.0483	
Std Dev	0.0006	
Rel Std Dev(%)	1.1945	


Operator's Signature

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001255
05/08/2026
Software: 8100.27

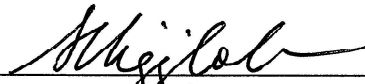
Test	g/210L	Time
Air Blank	0.000	12:43
Control Test	0.078	12:43
Air Blank	0.000	12:44
Control Test	0.077	12:44
Air Blank	0.000	12:45
Control Test	0.077	12:46
Air Blank	0.000	12:46
Control Test Stats		
Average	0.0773 ✓	
Std Dev	0.0006	
Rel Std Dev(%)	0.7466	


Operator's Signature

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001255
05/08/2026
Software: 8100.27

Test	g/210L	Time

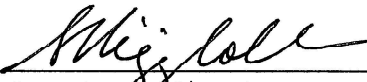
Air Blank	0.000	12:30
Control Test	0.198	12:31
Air Blank	0.000	12:31
Control Test	0.198	12:32
Air Blank	0.000	12:32
Control Test	0.197	12:33
Air Blank	0.000	12:34
Control Test Stats		
Average	0.1977	
Std Dev	0.0006	
Rel Std Dev(%)	0.2921	


Operator's Signature

ORANGE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001255
05/08/2026
Software: 8100.27

DGS

Test	g/210L	Time
Air Blank	0.000	13:12
Control Test	0.080	13:12
Air Blank	0.000	13:13
Control Test	0.080	13:13
Air Blank	0.000	13:14
Control Test	0.080	13:14
Air Blank	0.000	13:14
Control Test Stats		
Average	0.0800 ✓	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	


Operator Signature

instrument able to
read DGS pressure ✓
SWS

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: ORANGE COUNTY SO
Time of Inspection: 13:10

Date of Inspection: 05/18/2026

Serial Number: 80-001255
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202603D Exp: 03/04/2028	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:23825080A1 Exp: 10/05/2027
0.000	0.049	0.078	0.199	0.080
0.000	0.049	0.078	0.199	0.080
0.000	0.049	0.078	0.198	0.080
0.000	0.049	0.078	0.199	0.080
0.000	0.049	0.078	0.198	0.080
0.000	0.050	0.078	0.198	0.079
0.000	0.050	0.079	0.199	0.080
0.000	0.049	0.079	0.198	0.080
0.000	0.050	0.079	0.198	0.080
0.000	0.050	0.078	0.199	0.079

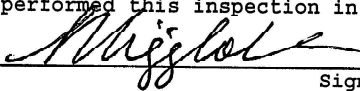
Standard Deviations	0.0005	0.0004	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.


LEANDRA HIGGINBOTHAM
Signature and Printed Name

05/18/2026
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001255, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001255</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>ORANGE COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>05/18/2026</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>13:10</u>	0.200 g/ 210 L	0.008
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

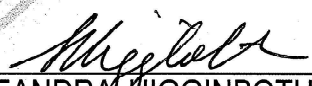
Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/18/2026

Date


LEANDRA HIGGINBOTHAM,
Department Inspector

FDLE/ATP Form 69 January 2026

Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

Page 1 of 1