



INSTRUMENT PROCESSING SHEET

Agency SEMINOLE COUNTY SOS/N 80-001271Florida Department of
Law EnforcementDate In 9/24/25DI Completion Date 10/2/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By <u>KTS</u>	Date <u>9/24/25</u>	Quality Checks	By <u>WKP</u>	Date <u>09/25/2025</u>	Flow Calibration	By _____	Date _____																																					
☒ Annual ☐ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: Intoxilyzer box came inside another padded FedEx box. KTS 9/24/25			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>132</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 103</u> 32 mm <u>0.148</u> (.139 - .169) 36 mm <u>0.164</u> (.156 - .190) 53 mm <u>0.238</u> (.228 - .278) 103 mm <u>0.500</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28421</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min - 17mm ☐ 15L/min - 53mm ☐ 30L/min - 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																							
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG510701 04/17/2027</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202406K 06/19/2026	0.080	MP5089	202406L 06/19/2026	0.200	MP5090	202406N 06/20/2026	0.080 DGS	N/A	AG510701 04/17/2027	<table border="1"><thead><tr><th colspan="2">Maintenance By _____ Date _____</th></tr><tr><td>☐ Battery Replacement</td><td></td></tr><tr><td>☐ Dry Gas Regulator Replacement</td><td></td></tr><tr><td>☐ Breath Tube Replacement</td><td></td></tr><tr><td>☐ Other _____</td><td></td></tr></thead></table>			Maintenance By _____ Date _____		☐ Battery Replacement		☐ Dry Gas Regulator Replacement		☐ Breath Tube Replacement		☐ Other _____													
Simulator	Serial #	Lot #/Exp																																											
0.050	MP5088	202406K 06/19/2026																																											
0.080	MP5089	202406L 06/19/2026																																											
0.200	MP5090	202406N 06/20/2026																																											
0.080 DGS	N/A	AG510701 04/17/2027																																											
Maintenance By _____ Date _____																																													
☐ Battery Replacement																																													
☐ Dry Gas Regulator Replacement																																													
☐ Breath Tube Replacement																																													
☐ Other _____																																													
Calibration Adjustment By <u>SLH</u>			Department Inspection By <u>SLH</u>																																										
Barometric Pressure Gauge <u>1008</u> ID # <u>28662</u>			Barometric Pressure ID# <u>28421</u> Gauge <u>1017</u> Instrument <u>1016</u> Mouth Alcohol Solution Lot # <u>2025-C</u> Acetone Stock Solution Lot # <u>2024-B</u>																																										
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td>MP6294</td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td>MP6295</td><td>24330</td><td>09/10/2026</td></tr><tr><td>0.100</td><td>MP6296</td><td>24390</td><td>10/29/2026</td></tr><tr><td>0.200</td><td>MP6297</td><td>24080</td><td>02/13/2026</td></tr><tr><td>0.300</td><td>MP6298</td><td>24430</td><td>12/10/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>28424080A3</td><td>11/05/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000	MP6294	N/A	N/A	0.040	MP6295	24330	09/10/2026	0.100	MP6296	24390	10/29/2026	0.200	MP6297	24080	02/13/2026	0.300	MP6298	24430	12/10/2026	0.080 DGS	N/A	28424080A3	11/05/2026	<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6294</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP5088</td></tr><tr><td>0.080</td><td>MP5089</td></tr><tr><td>0.200</td><td>MP5090</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP6294	Interferent	MP6290	0.050	MP5088	0.080	MP5089	0.200	MP5090
Simulator	Serial #	Lot #	Expiration																																										
0.000	MP6294	N/A	N/A																																										
0.040	MP6295	24330	09/10/2026																																										
0.100	MP6296	24390	10/29/2026																																										
0.200	MP6297	24080	02/13/2026																																										
0.300	MP6298	24430	12/10/2026																																										
0.080 DGS	N/A	28424080A3	11/05/2026																																										
Simulator	Serial Number																																												
0.000	MP6294																																												
Interferent	MP6290																																												
0.050	MP5088																																												
0.080	MP5089																																												
0.200	MP5090																																												
☒ Post Calibration Adjustment Stability Checks			Attachments																																										
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td>MP6291</td><td>202406K</td><td>06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L</td><td>06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N</td><td>06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>** 28424080A3</td><td>10/22/2026</td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050	MP6291	202406K	06/19/2026	0.080	MP5089	202406L	06/19/2026	0.200	MP5090	202406N	06/20/2026	0.080 DGS	N/A	** 28424080A3	10/22/2026	<table border="1"><tbody><tr><td>☒ Form 41</td><td>☒ Post-Stability Checks</td></tr><tr><td>☒ Stability Checks</td><td>☐ Flow Calibration</td></tr><tr><td>☒ Calibration Certificate</td><td>☐ Form 40</td></tr><tr><td>☒ Calibration Adjustment</td><td>☒ Other <u>added stabilities</u></td></tr></tbody></table>			☒ Form 41	☒ Post-Stability Checks	☒ Stability Checks	☐ Flow Calibration	☒ Calibration Certificate	☐ Form 40	☒ Calibration Adjustment	☒ Other <u>added stabilities</u>												
Simulator	Serial #	Lot #	Expiration																																										
0.050	MP6291	202406K	06/19/2026																																										
0.080	MP5089	202406L	06/19/2026																																										
0.200	MP5090	202406N	06/20/2026																																										
0.080 DGS	N/A	** 28424080A3	10/22/2026																																										
☒ Form 41	☒ Post-Stability Checks																																												
☒ Stability Checks	☐ Flow Calibration																																												
☒ Calibration Certificate	☐ Form 40																																												
☒ Calibration Adjustment	☒ Other <u>added stabilities</u>																																												
Notes/Suggested Service: The average of the 0.08 wet and the average of the 0.08 DGS are more than 0.003 apart. Post stabilities performed 2 separate days. SLH 9/30/25 *Tech Correction: Corrected missing DGS first stability page. WKP 10/7/2025 ** Tech Review- Corrected lot # for DGS in post-stabilities- correct lot # is AG429602. Also populated post-calibration stability box. SLH 10/7/2025			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use Taylor Gutschow Digitally signed by Taylor Gutschow Date: 2025.10.07 16:07:37 -04'00' Platt Digitally signed by Shayla Platt Date: 2025.10.15 19:10:07 -04'00' Tech Review / Date _____ Admin Review / Date _____																																										

Stability Checks

0.05g/210L	0.08g/210L	0.20g/210L	DGS 0.08g/210L																																																																																																																																																
0.047 to 0.053	0.077 to 0.083	0.194 to 0.206	0.077 to 0.083																																																																																																																																																
✓	✓	✓	✓																																																																																																																																																
Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis																																																																																																																																																
✓	✓	✓	✓																																																																																																																																																
SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/25/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>11:38</td></tr> <tr><td>Control Test</td><td>0.050</td><td>11:39</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:40</td></tr> <tr><td>Control Test</td><td>0.048</td><td>11:40</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:41</td></tr> <tr><td>Control Test</td><td>0.049</td><td>11:42</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:42</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0490</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>2.0408</td><td></td></tr> </tbody> </table> Operator's Signature: <i>menpi</i>	Test	g/210L	Time	Air Blank	0.000	11:38	Control Test	0.050	11:39	Air Blank	0.000	11:40	Control Test	0.048	11:40	Air Blank	0.000	11:41	Control Test	0.049	11:42	Air Blank	0.000	11:42	Control Test Stats			Average	0.0490		Std Dev	0.0010		Rel Std Dev(%)	2.0408		SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/25/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>11:48</td></tr> <tr><td>Control Test</td><td>0.079</td><td>11:49</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:50</td></tr> <tr><td>Control Test</td><td>0.079</td><td>11:50</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:51</td></tr> <tr><td>Control Test</td><td>0.078</td><td>11:52</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:52</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0787</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7339</td><td></td></tr> </tbody> </table> Operator's Signature: <i>menpi</i>	Test	g/210L	Time	Air Blank	0.000	11:48	Control Test	0.079	11:49	Air Blank	0.000	11:50	Control Test	0.079	11:50	Air Blank	0.000	11:51	Control Test	0.078	11:52	Air Blank	0.000	11:52	Control Test Stats			Average	0.0787		Std Dev	0.0006		Rel Std Dev(%)	0.7339		SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/25/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>11:54</td></tr> <tr><td>Control Test</td><td>0.201</td><td>11:55</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:55</td></tr> <tr><td>Control Test</td><td>0.199</td><td>11:56</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:57</td></tr> <tr><td>Control Test</td><td>0.200</td><td>11:57</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:58</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.2000</td><td></td></tr> <tr><td>Std Dev</td><td>0.0010</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.5000</td><td></td></tr> </tbody> </table> Operator's Signature: <i>menpi</i>	Test	g/210L	Time	Air Blank	0.000	11:54	Control Test	0.201	11:55	Air Blank	0.000	11:55	Control Test	0.199	11:56	Air Blank	0.000	11:57	Control Test	0.200	11:57	Air Blank	0.000	11:58	Control Test Stats			Average	0.2000		Std Dev	0.0010		Rel Std Dev(%)	0.5000		SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/25/2025 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:01</td></tr> <tr><td>Control Test</td><td>0.082</td><td>14:01</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:01</td></tr> <tr><td>Control Test</td><td>0.082</td><td>14:02</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:02</td></tr> <tr><td>Control Test</td><td>0.082</td><td>14:02</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:02</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0820</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> Operator's Signature: <i>menpi</i>	Test	g/210L	Time	Air Blank	0.000	14:01	Control Test	0.082	14:01	Air Blank	0.000	14:01	Control Test	0.082	14:02	Air Blank	0.000	14:02	Control Test	0.082	14:02	Air Blank	0.000	14:02	Control Test Stats			Average	0.0820		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:38																																																																																																																																																	
Control Test	0.050	11:39																																																																																																																																																	
Air Blank	0.000	11:40																																																																																																																																																	
Control Test	0.048	11:40																																																																																																																																																	
Air Blank	0.000	11:41																																																																																																																																																	
Control Test	0.049	11:42																																																																																																																																																	
Air Blank	0.000	11:42																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0490																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	2.0408																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:48																																																																																																																																																	
Control Test	0.079	11:49																																																																																																																																																	
Air Blank	0.000	11:50																																																																																																																																																	
Control Test	0.079	11:50																																																																																																																																																	
Air Blank	0.000	11:51																																																																																																																																																	
Control Test	0.078	11:52																																																																																																																																																	
Air Blank	0.000	11:52																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0787																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7339																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:54																																																																																																																																																	
Control Test	0.201	11:55																																																																																																																																																	
Air Blank	0.000	11:55																																																																																																																																																	
Control Test	0.199	11:56																																																																																																																																																	
Air Blank	0.000	11:57																																																																																																																																																	
Control Test	0.200	11:57																																																																																																																																																	
Air Blank	0.000	11:58																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.2000																																																																																																																																																		
Std Dev	0.0010																																																																																																																																																		
Rel Std Dev(%)	0.5000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:01																																																																																																																																																	
Control Test	0.082	14:01																																																																																																																																																	
Air Blank	0.000	14:01																																																																																																																																																	
Control Test	0.082	14:02																																																																																																																																																	
Air Blank	0.000	14:02																																																																																																																																																	
Control Test	0.082	14:02																																																																																																																																																	
Air Blank	0.000	14:02																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0820																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		

SEMINOLE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 60-001271
09/25/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	11:33
Control Test	INT*	11:34
Air Blank	0.000	11:35
Control Test	0.049	11:36
Air Blank	0.000	11:36
Control Test	0.049	11:37
Air Blank	0.000	11:37
Control Test Stats		
Average	0.0497	
Std Dev	0.0012	
Rel Std Dev(%)	2.3249	

*Interferent Detect


Operator's Signature

0.05 wet stability

Inadvertently unplugged simulator during the first simulator, then ^{run} plugged it back in for the second and third run. Repeated analysis post root cause analysis.

WKP 09/25/2025

SEMINOLE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001271
09/25/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	12:10
Control Test	0.082	12:10
Air Blank	0.000	12:10
Control Test	0.082	12:11
Air Blank	0.000	12:11
Control Test	INT*	12:12
Air Blank	0.000	12:12
Control Test Stats		
Average	0.0820	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

*Interferent Detect

menpi
Operator's Signature

0.08g/210L
DBS

INT* exception

message prompted.

20 stabilities of
the DBS were
performed to
determine if *INT
detect could be
reproduced.

wkp 09/25/2025

SEMINOLE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001271
09/25/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:19
Control Test	0.049	13:20
Air Blank	0.000	13:20
Control Test	0.048	13:21
Air Blank	0.000	13:22
Control Test	0.049	13:22
Air Blank	0.000	13:23
Control Test	0.048	13:23
Air Blank	0.000	13:24
Control Test	0.048	13:25
Air Blank	0.000	13:25
Control Test	0.048	13:26
Air Blank	0.000	13:27
Control Test	0.048	13:27
Air Blank	0.000	13:28
Control Test	0.047	13:29
Air Blank	0.000	13:29
Control Test	0.047	13:30
Air Blank	0.000	13:30
Control Test	0.047	13:31
Air Blank	0.000	13:32
Control Test	0.048	13:32
Air Blank	0.000	13:33
Control Test	0.047	13:34
Air Blank	0.000	13:34
Control Test	0.047	13:35
Air Blank	0.000	13:35
Control Test	0.047	13:36
Air Blank	0.000	13:37
Control Test	0.048	13:37
Air Blank	0.000	13:38
Control Test	0.048	13:39
Air Blank	0.000	13:39
Control Test	0.047	13:40
Air Blank	0.000	13:40
Control Test	0.047	13:41
Air Blank	0.000	13:42
Control Test	0.048	13:42
Air Blank	0.000	13:43
Control Test	0.047	13:44
Air Blank	0.000	13:44
Control Test Stats		
Average	0.0476	
Std Dev	0.0007	
Rel Std Dev(%)	1.4078	

Operator's Signature

20 stabilities
Performed to
determine if

INT* could
be reproduced
for both
0.05 g/210L
and 0.08 g/210L
DGS. Stabilities
were reperfomed
for the corresponding
concentrations.

09/25/2025

SEMINOLE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001271
09/25/2025
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	13:01
Control Test	0.082	13:02
Air Blank	0.000	13:02
Control Test	0.082	13:02
Air Blank	0.000	13:03
Control Test	0.082	13:03
Air Blank	0.000	13:04
Control Test	0.082	13:04
Air Blank	0.000	13:05
Control Test	0.082	13:05
Air Blank	0.000	13:05
Control Test	0.082	13:06
Air Blank	0.000	13:06
Control Test	0.082	13:07
Air Blank	0.000	13:07
Control Test	0.082	13:07
Air Blank	0.000	13:08
Control Test	0.082	13:08
Air Blank	0.000	13:09
Control Test	0.083	13:09
Air Blank	0.000	13:10
Control Test	0.083	13:10
Air Blank	0.000	13:10
Control Test	0.082	13:11
Air Blank	0.000	13:11
Control Test	0.082	13:11
Air Blank	0.000	13:12
Control Test	0.081	13:12
Air Blank	0.000	13:13
Control Test	0.082	13:13
Air Blank	0.000	13:14
Control Test	0.082	13:14
Air Blank	0.000	13:14
Control Test	0.082	13:15
Air Blank	0.000	13:15
Control Test	0.082	13:16
Air Blank	0.000	13:16
Control Test	0.082	13:16
Air Blank	0.000	13:17
Control Test	0.081	13:17
Air Blank	0.000	13:18
Control Test Stats		
Average	0.0820	
Std Dev	0.0005	
Rel Std Dev(%)	0.5595	

Operator's Signature

SEMINOLE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001271
09/30/2025 13:42:48

80-001271
SUN 9/30/25

Auto Calibration

pg 1 of 2

<<<< 3um >>>> <<<< 9um >>>>

Solution = 0.000 g/210L or 0.0000 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	0.1350	(-0.0160)	0.1960	(-0.0100)
Sample #2	0.0430	(0.1020)	0.1840	(-0.0070)
Sample #3	0.1350	(0.1030)	0.1760	(0.0130)
Sample #4	0.0740	(0.0930)	0.1700	(0.0190)
Avg % Abs	0.0840	(0.0993)	0.1767	(0.0083)
STD DEV	0.0468	(0.0055)	0.0070	(0.0136)
REL STD DEV	55.724	(5.545)	3.976	(163.365)

Solution = 0.040 g/210L or 0.1905 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	0.8350	(-0.0170)	1.5450	(-0.0040)
Sample #2	0.8450	(0.0400)	1.5150	(0.0420)
Sample #3	0.8490	(0.0540)	1.5230	(0.0490)
Sample #4	0.8250	(0.0600)	1.5140	(0.0520)
Avg % Abs	0.8397	(0.0513)	1.5173	(0.0477)
STD DEV	0.0129	(0.0103)	0.0049	(0.0051)
REL STD DEV	1.531	(19.993)	0.325	(10.766)

Solution = 0.100 g/210L or 0.4762 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	1.9070	(-0.0130)	3.5890	(0.0040)
Sample #2	1.8750	(0.0250)	3.5340	(0.0460)
Sample #3	1.8660	(0.0620)	3.5080	(0.0700)
Sample #4	1.8580	(0.0600)	3.4950	(0.0800)
Avg % Abs	1.8663	(0.0490)	3.5123	(0.0653)
STD DEV	0.0085	(0.0208)	0.0199	(0.0175)
REL STD DEV	0.456	(42.467)	0.565	(26.746)

Solution = 0.200 g/210L or 0.9524 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	3.6820	(-0.0300)	6.8040	(-0.0150)
Sample #2	3.6000	(0.0680)	6.7000	(0.0810)
Sample #3	3.5760	(0.1050)	6.6820	(0.1140)
Sample #4	3.5720	(0.1170)	6.6640	(0.1270)
Avg % Abs	3.5827	(0.0967)	6.6820	(0.1073)
STD DEV	0.0151	(0.0255)	0.0180	(0.0237)
REL STD DEV	0.423	(26.422)	0.269	(22.093)

Solution = 0.300 g/210L or 1.4286 mg/l, Samples = 4, Discarded = 1

Sample	% Abs	(% Abs Ref)	% Abs	(% Abs Ref)
Sample #1	5.3470	(-0.0060)	9.8270	(0.0020)
Sample #2	5.2560	(0.1240)	9.7180	(0.1640)
Sample #3	5.2830	(0.1130)	9.7030	(0.1670)
Sample #4	5.2530	(0.1430)	9.6780	(0.1920)
Avg % Abs	5.2640	(0.1267)	9.6997	(0.1743)
STD DEV	0.0165	(0.0152)	0.0202	(0.0154)
REL STD DEV	0.314	(11.982)	0.208	(8.818)

SEMINOLE COUNTY SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-001271
09/30/2025 13:42:48

80-001271
Sub 9/30/25

Auto Calibration

pg 2 of 2

<<<< 3um >>>>

Zero Order Coef -261.48
First Order Coef 2638.67
Second Order Coef 24.10

<<<< 9um >>>>

-235.26
1382.35
11.80

Act Fit Residual
(g/210L) (g/210L) (g/210L)
0.000 -0.001 0.0008
0.040 0.041 -0.0014
0.100 0.100 0.0003
0.200 0.200 0.0005
0.300 0.300 -0.0002

Act Fit Residual
(g/210L) (g/210L) (g/210L)
0.000 0.000 -0.0002
0.040 0.040 0.0003
0.100 0.100 -0.0001
0.200 0.200 -0.0001
0.300 0.300 0.0001

<<<< 3um >>>> <<<< 9um >>>>

Solution = 0.080 g/210L or 0.3810 mg/l, Samples = 4, Discarded = 1
Sample
Sample #1 3329.00 3306.00
Sample #2 3348.00 3322.00
Sample #3 3247.00 3335.00
Sample #4 3217.00 3311.00
Avg 3270.6667 3322.6667
STD DEV 68.6319 12.0139
REL STD DEV 2.098 0.362
H2O adjust (mg/l*10k) 539 487

Barometric Pressure = 1008

*****CALIBRATION SUCCESSFUL*****

Post-Calibration Adjustment Stability Checks

80-001271
*9/30/2025
S44

* 0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	* DGS 0.080 g/210L 0.077 to 0.083 g/210L ±0.003 g/210L of Wet
SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/30/2025 Software: 8100.27 SN 80-001271 Test g/210L Time Air Blank 0.010 14:56 Control Test 0.050 14:57 Air Blank 0.010 14:57 Control Test 0.049 14:58 Air Blank 0.010 14:59 Control Test 0.050 14:59 Air Blank 0.000 15:00 Control Test Stats Average 0.0497 Std Dev 0.0006 Rel Std Dev(%) 1.1625	SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 10/02/2025 Software: 8100.27 SN 80-001271 Test g/210L Time Air Blank 0.000 12:17 Control Test 0.079 12:17 Air Blank 0.000 12:18 Control Test 0.080 12:18 Air Blank 0.000 12:19 Control Test 0.079 12:20 Air Blank 0.000 12:20 Control Test Stats Average 0.0793 Std Dev 0.0006 Rel Std Dev(%) 0.7277	SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 10/02/2025 Software: 8100.27 SN 80-001271 Test g/210L Time Air Blank 0.000 11:58 Control Test 0.203 11:59 Air Blank 0.100 11:59 Control Test 0.201 12:00 Air Blank 0.000 12:01 Control Test 0.200 12:01 Air Blank 0.000 12:02 Control Test Stats Average 0.2113 Std Dev 0.0015 Rel Std Dev(%) 0.7587	SEMINOLE COUNTY SO Intoxilyzer - Alcohol Analyzer Model 8000 09/30/2025 Software: 8100.27 SN 80-001271 Test g/210L Time Air Blank 0.000 14:48 Control Test 0.079 14:48 Air Blank 0.000 14:49 Control Test 0.079 14:49 Air Blank 0.000 14:50 Control Test 0.079 14:50 Air Blank 0.000 14:50 Control Test Stats Average 0.0790 Std Dev 0.0000 Rel Std Dev(%) 0.0000
Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis	Performed Root Case Analysis
Operator's Signature	Operator's Signature	Operator's Signature	Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: SEMINOLE COUNTY SO
Time of Inspection: 16:23

Date of Inspection: 10/02/2025

Serial Number: 80-001271
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG510701 Exp: 04/17/2027
0.000	0.049	0.080	0.201	0.078
0.000	0.050	0.080	0.201	0.078
0.000	0.050	0.080	0.200	0.078
0.000	0.050	0.080	0.200	0.078
0.000	0.050	0.080	0.201	0.078
0.000	0.050	0.080	0.200	0.077
0.000	0.050	0.080	0.201	0.077
0.000	0.050	0.081	0.201	0.078
0.000	0.050	0.080	0.201	0.078
0.000	0.049	0.081	0.201	0.078
Standard Deviations	0.0004	0.0004	0.0004	0.0004

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5



Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

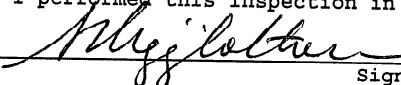
LEANDRA HIGGINBOTHAM
Signature and Printed Name

10/02/2025
Date

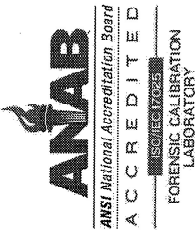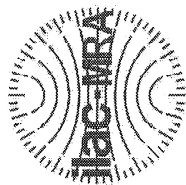

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001271, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001271</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>SEMINOLE COUNTY SO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>10/02/2025</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>16:23</u>	0.200 g/ 210 L	0.007
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 October 2024

Issuing Authority: Alcohol Testing Program

10/02/2025

Date

Mejlat

LEANDRA HIGGINBOTHAM,

Department Inspector

Service • Integrity • Respect • Quality