



INSTRUMENT PROCESSING SHEET

Agency Miami Dade SOS/N 80-000884Florida Department of
Law EnforcementDate In 9/8/2025DI Completion Date 09/16/2025☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake	By SLH	Date 9/8/2025	Quality Checks	By WKP	Date 9/9/2025	Flow Calibration	By	Date																														
☒ Annual ☒ Registration ☐ Return from CMI / EE Visual Inspection: ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/ Accessories: ☐ Power cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes: <u>Unattached foot-included.</u> <u>Note to change name to MDSO.</u>			☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value <u>147</u> ☒ Flow Verification (L/s) Flow Column # <u>ATP 102</u> 32 mm <u>0.144</u> (.139 - .169) 36 mm <u>0.160</u> (.156 - .190) 53 mm <u>0.230</u> (.228 - .278) 103 mm <u>0.496</u> (.447 - .547) ☒ Barometric Pressure Check Gauge ID # <u>28662</u> ☒ Stability Checks			Flow Column # _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value _____ ☐ Post Calibration Verification (L/s) Flow Column # _____ 32 mm _____ (.139 - .169) 36 mm _____ (.156 - .190) 53 mm _____ (.228 - .278) 103 mm _____ (.447 - .547)																																
			<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202406K 06/19/2026</td></tr><tr><td>0.080</td><td>MP5089</td><td>202406L 06/19/2026</td></tr><tr><td>0.200</td><td>MP5090</td><td>202406N 06/20/2026</td></tr><tr><td>0.080 DGS</td><td>N/A</td><td>AG510701 04/17/2027</td></tr></tbody></table>			Simulator	Serial #	Lot #/Exp	0.050	MP5088	202406K 06/19/2026	0.080	MP5089	202406L 06/19/2026	0.200	MP5090	202406N 06/20/2026	0.080 DGS	N/A	AG510701 04/17/2027	<table border="1"><thead><tr><th colspan="2">Maintenance By SLH</th><th>Date 9/10/25</th></tr></thead><tbody><tr><td colspan="2">☐ Battery Replacement</td><td></td></tr><tr><td colspan="2">☐ Dry Gas Regulator Replacement</td><td></td></tr><tr><td colspan="2">☐ Breath Tube Replacement</td><td></td></tr><tr><td colspan="2">☒ Other <u>Attached foot</u></td><td></td></tr></tbody></table>			Maintenance By SLH		Date 9/10/25	☐ Battery Replacement			☐ Dry Gas Regulator Replacement			☐ Breath Tube Replacement			☒ Other <u>Attached foot</u>		
Simulator	Serial #	Lot #/Exp																																				
0.050	MP5088	202406K 06/19/2026																																				
0.080	MP5089	202406L 06/19/2026																																				
0.200	MP5090	202406N 06/20/2026																																				
0.080 DGS	N/A	AG510701 04/17/2027																																				
Maintenance By SLH		Date 9/10/25																																				
☐ Battery Replacement																																						
☐ Dry Gas Regulator Replacement																																						
☐ Breath Tube Replacement																																						
☒ Other <u>Attached foot</u>																																						
Calibration Adjustment			Department Inspection																																			
By _____			By <u>SLH</u>																																			
Barometric Pressure Gauge _____ ID # _____			Barometric Pressure ID# <u>28421</u>																																			
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge <u>1016</u> Instrument <u>1015</u> Mouth Alcohol Solution Lot # <u>2025-A</u> Acetone Stock Solution Lot # <u>2024-B</u>							
Simulator	Serial #	Lot #	Expiration																																			
0.000		N/A	N/A																																			
0.040																																						
0.100																																						
0.200																																						
0.300																																						
0.080 DGS	N/A																																					
☐ Post Calibration Adjustment Stability Checks			<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP5087</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>			Simulator	Serial Number	0.000	MP6289	Interferent	MP5087	0.050	MP6291	0.080	MP6292	0.200	MP6293																					
Simulator	Serial Number																																					
0.000	MP6289																																					
Interferent	MP5087																																					
0.050	MP6291																																					
0.080	MP6292																																					
0.200	MP6293																																					
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>			Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			Attachments ☒ Form 41 ☐ Post-Stability Checks ☒ Stability Checks ☐ Flow Calibration ☒ Calibration Certificate ☐ Form 40 ☐ Calibration Adjustment ☒ Other <u>Note&Diagnostic</u>															
Simulator	Serial #	Lot #	Expiration																																			
0.050																																						
0.080																																						
0.200																																						
0.080 DGS	N/A																																					
Notes/Suggested Service: _____			☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																			
			Shayla Platt Digitally signed by Shayla Platt Date: 2025.10.02 08:03:17 -04'00'																																			
			Tech Review / Date _____ Admin Review / Date _____																																			

Taylor,

80-000884 —

Standard Annual.

Needs to be changed

to

MDSO.

Thx,

— R. Closius

Rec'd 9/8/25
SWA

Stability Checks

Noted-Agency name not updated for stabilities.
SLH 9/18/2025

0.05g/210L			0.08g/210L			0.20g/210L			DGS 0.08g/210L		
0.047 to 0.053			0.077 to 0.083			0.194 to 0.206			0.077 to 0.083		
Performed Root Case Analysis			Performed Root Case Analysis			Performed Root Case Analysis			Performed Root Case Analysis		
MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model: 8000 SN 80-000664 09/09/2025 Software: 8100.27			MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model: 8000 SN 80-000664 09/09/2025 Software: 8100.27			MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model: 8000 SN 80-000664 09/09/2025 Software: 8100.27			MIAMI DADE PD Intoxilyzer - Alcohol Analyzer Model: 8000 SN 80-000664 09/09/2025 Software: 8100.27		
Test	g/210L	Time	Test	g/210L	Time	Test	g/210L	Time	Test	g/210L	Time
Air Blank	0.000	08:40	Air Blank	0.000	08:45	Air Blank	0.000	08:51	Air Blank	0.000	08:56
Control Test	0.049	08:41	Control Test	0.080	08:46	Control Test	0.199	08:51	Control Test	0.080	08:56
Air Blank	0.000	08:42	Air Blank	0.000	08:47	Air Blank	0.000	08:52	Air Blank	0.000	08:58
Control Test	0.049	08:42	Control Test	0.079	08:47	Control Test	0.198	08:53	Control Test	0.081	08:59
Air Blank	0.000	08:43	Air Blank	0.000	08:48	Air Blank	0.000	08:53	Air Blank	0.000	08:59
Control Test	0.049	08:44	Control Test	0.079	08:49	Control Test	0.198	08:54	Control Test	0.082	09:00
Air Blank	0.000	08:44	Air Blank	0.000	08:49	Air Blank	0.000	08:54	Air Blank	0.000	09:00
Control Test Stats			Control Test Stats			Control Test Stats			Control Test Stats		
Average	0.0490		Average	0.0793		Average	0.1960		Average	0.0813	
Std Dev	0.0000		Std Dev	0.0006		Std Dev	0.0000		Std Dev	0.0006	
Rel Std Dev(%)	0.0000		Rel Std Dev(%)	0.7277		Rel Std Dev(%)	0.0000		Rel Std Dev(%)	0.7099	
Operator's Signature			Operator's Signature			Operator's Signature			Operator's Signature		

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

9/16/25 SUT
Agency: MIAMI DADE PD **SO**
Time of Inspection: 17:08

Date of Inspection: 09/16/2025

Serial Number: 80-000884
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202406K Exp: 06/19/2026	0.08g/210L Test (g/210L) Lot#:202406L Exp: 06/19/2026	0.20g/210L Test (g/210L) Lot#:202406N Exp: 06/20/2026	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:AG510701 Exp: 04/17/2027
0.000	0.050	0.079	0.199	0.081
0.000	0.050	0.079	0.199	0.081
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.079	0.199	0.080
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.079	0.199	0.079
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.079	0.200	0.079
0.000	0.050	0.080	0.199	0.079
0.000	0.050	0.079	0.199	0.079


Standard Deviations	0.0000	0.0004	0.0003	0.0008
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0003 Number of Simulators Used: 5


Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

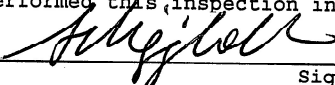
LEANDRA HIGGINBOTHAM
Signature and Printed Name

09/16/2025
Date

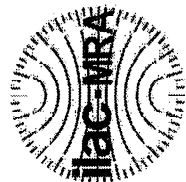

ANAB
ANSI National Accreditation Board
ACCREDITED
ISO/IEC 17025
FORENSIC CALIBRATION
LABORATORY

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-000884, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	80-000884	UNCERTAINTY* $\pm$
Owning Agency:	MIAMI DADE PD	0.050 g/ 210 L
Calibration Date:	09/16/2025	0.080 g/ 210 L
Calibration Time:	17:08	0.200 g/ 210 L
		0.080 g/ 210 L Dry Gas Control
		0.004
		0.004
		0.007
		0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACCS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards. This document shall not be reproduced except in full without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

FDLE/ATP Form 69 October 2024

Issuing Authority: Alcohol Testing Program

Date

09/16/2025

LEANDRA HIGGINBOTHAM,

Department Inspector

Service • Integrity • Respect • Quality

MIAMI DADE SO
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-000884
09/16/2025
Software: 8100.27

: DIAGNOSTICS :

Voltage/Current Test	OK
RAM Test	OK
EEPROM Checksum Test	OK
Real Time Clock Test	OK
DSP Test	OK
Analytical Stability Test	OK
Internal Printer Test	OK
Modem Test	OK
Temperature Regulation Test	OK

* Shows
- named changed - Sub
post inspection
and calibration
certificate generation
9/16/25