

Alcohol Testing Program - Instrument Processing Sheet

Agency: LAKELAND PD

Instrument Serial Number: 80-003945

Date In: 7/7/2026

DI Completion Date: 7/16/2026

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By: <u>WKP</u> Date: <u>7/7/2026</u>				Quality Checks By: <u>KTS</u> Date: <u>7/16/2026</u>				Flow Adjustment By: _____			
☒ Annual ☐ Dropped Off ☐ Registration ☐ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☒ Power Cord ☐ Printer Cable ☐ Static Bag ☐ 12V DC Cable Notes: Instrument arrived with two back dry gas shelf feet missing. WKP 7/7/26				☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>204</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP103</u> 32 mm <u>0.152</u> (.139-.169) 36 mm <u>0.164</u> (.156-.190) 53 mm <u>0.234</u> (.228-.278) 103 mm <u>0.492</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1012</u> Instrument: <u>1005</u> ☒ Stability Checks				Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)			
				Simulator		Serial #		Lot#/Exp		Maintenance By: <u>KTS</u> Date: <u>7/16/2026</u>	
				0.050		MP5088		202603A 3/3/2028		☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Tank Sensor Tare ☐ Breath Tube Replacement ☒ Other: Added electrical tape to breath hose.	
				0.080		MP5089		202603B 3/3/2028			
				0.200		MP5090		202603D 3/4/2028			
0.080 DGS		N/A		23825080A1 10/5/2027							
Optical Bench Adjustment By: _____				Department Inspection By: <u>KTS</u>							
Barometric Pressure Gauge: _____ ID#: _____				Barometric Pressure ID#: <u>28427</u>							
Simulator		Serial #		Lot #		Expiration		Gauge: <u>1014</u> Instrument: <u>1007</u>			
0.000				N/A		N/A		Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>9/25/2027</u>			
0.040								Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>9/22/2027</u>			
0.100								Simulator			
0.200								Serial Number			
0.300								0.000			
0.080 DGS		N/A						MP5086			
☐ Post Optical Bench Adjustment Stability Checks								Interferent			
Simulator		Serial #		Lot #		Expiration		MP5087			
0.050								0.050			
0.080								MP5088			
0.200								0.080			
0.080 DGS		N/A						MP5089			
Gauge ID #: _____				Gauge: _____ Instrument: _____				0.200			
								MP5090			
Notes/Suggested Service:				Attachments							
				☒ Form 41				☐ Post-Stability Checks			
				☒ Stability Checks				☐ Flow Adjustment			
				☒ Calibration Certificate				☐ Form 40			
				☐ Optical Bench Adjustment				☐ Other:			
				☒ Instrument Complies with Chapter 11D-8, FAC							
				☐ Instrument Does Not Comply with Chapter 11D-8, FAC							
				☒ Return to/Place into Evidentiary Use							
				☐ Remain Out of Evidentiary Use							
				☒ Conduct an Agency Inspection Before Evidentiary Use							
				LeAndra Higginbotham		Digitally signed by LeAndra Higginbotham Date: 2026.08.03 12:35:06 -04'00'		Wen-Chi Pierson			
				Tech Review		Admin Review		Digitally signed by Wen-Chi Pierson Date: 2026.08.04 08:53:29 -04'00'			

Stability Checks

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L ≤0.003 g/210L of Wet																																																																																																																																																
LAKELAND PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-003945 07/16/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>09:02</td></tr> <tr><td>Control Test</td><td>0.052</td><td>08:03</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>08:07</td></tr> <tr><td>Control Test</td><td>0.052</td><td>08:14</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>08:05</td></tr> <tr><td>Control Test</td><td>0.052</td><td>08:05</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>08:18</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0520</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> <i>Katz</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	09:02	Control Test	0.052	08:03	Air Blank	0.000	08:07	Control Test	0.052	08:14	Air Blank	0.000	08:05	Control Test	0.052	08:05	Air Blank	0.000	08:18	Control Test Stats			Average	0.0520		Std Dev	0.0000		Rel Std Dev(%)	0.0000		LAKELAND PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-003945 07/16/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>07:57</td></tr> <tr><td>Control Test</td><td>0.079</td><td>07:58</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:58</td></tr> <tr><td>Control Test</td><td>0.080</td><td>07:59</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>08:00</td></tr> <tr><td>Control Test</td><td>0.080</td><td>08:00</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>08:01</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0797</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7247</td><td></td></tr> </tbody> </table> <i>Katz</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	07:57	Control Test	0.079	07:58	Air Blank	0.000	07:58	Control Test	0.080	07:59	Air Blank	0.000	08:00	Control Test	0.080	08:00	Air Blank	0.000	08:01	Control Test Stats			Average	0.0797		Std Dev	0.0006		Rel Std Dev(%)	0.7247		LAKELAND PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-003945 07/16/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>07:50</td></tr> <tr><td>Control Test</td><td>0.199</td><td>07:51</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:51</td></tr> <tr><td>Control Test</td><td>0.198</td><td>07:52</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:52</td></tr> <tr><td>Control Test</td><td>0.198</td><td>07:53</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:54</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1983</td><td></td></tr> <tr><td>Std Dev</td><td>0.0008</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2911</td><td></td></tr> </tbody> </table> <i>Katz</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	07:50	Control Test	0.199	07:51	Air Blank	0.000	07:51	Control Test	0.198	07:52	Air Blank	0.000	07:52	Control Test	0.198	07:53	Air Blank	0.000	07:54	Control Test Stats			Average	0.1983		Std Dev	0.0008		Rel Std Dev(%)	0.2911		LAKELAND PD Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-003945 07/16/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>07:45</td></tr> <tr><td>Control Test</td><td>0.082</td><td>07:46</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:46</td></tr> <tr><td>Control Test</td><td>0.082</td><td>07:47</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:47</td></tr> <tr><td>Control Test</td><td>0.082</td><td>07:47</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>07:48</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0820</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> <i>Katz</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	07:45	Control Test	0.082	07:46	Air Blank	0.000	07:46	Control Test	0.082	07:47	Air Blank	0.000	07:47	Control Test	0.082	07:47	Air Blank	0.000	07:48	Control Test Stats			Average	0.0820		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:02																																																																																																																																																	
Control Test	0.052	08:03																																																																																																																																																	
Air Blank	0.000	08:07																																																																																																																																																	
Control Test	0.052	08:14																																																																																																																																																	
Air Blank	0.000	08:05																																																																																																																																																	
Control Test	0.052	08:05																																																																																																																																																	
Air Blank	0.000	08:18																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0520																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	07:57																																																																																																																																																	
Control Test	0.079	07:58																																																																																																																																																	
Air Blank	0.000	07:58																																																																																																																																																	
Control Test	0.080	07:59																																																																																																																																																	
Air Blank	0.000	08:00																																																																																																																																																	
Control Test	0.080	08:00																																																																																																																																																	
Air Blank	0.000	08:01																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0797																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7247																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	07:50																																																																																																																																																	
Control Test	0.199	07:51																																																																																																																																																	
Air Blank	0.000	07:51																																																																																																																																																	
Control Test	0.198	07:52																																																																																																																																																	
Air Blank	0.000	07:52																																																																																																																																																	
Control Test	0.198	07:53																																																																																																																																																	
Air Blank	0.000	07:54																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1983																																																																																																																																																		
Std Dev	0.0008																																																																																																																																																		
Rel Std Dev(%)	0.2911																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	07:45																																																																																																																																																	
Control Test	0.082	07:46																																																																																																																																																	
Air Blank	0.000	07:46																																																																																																																																																	
Control Test	0.082	07:47																																																																																																																																																	
Air Blank	0.000	07:47																																																																																																																																																	
Control Test	0.082	07:47																																																																																																																																																	
Air Blank	0.000	07:48																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0820																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: LAKELAND PD

Time of Inspection: 11:22

Date of Inspection: 07/16/2026

Serial Number: 80-003945

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202603D Exp: 03/04/2028	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:23825080A1 Exp: 10/05/2027
0.000	0.051	0.078	0.198	0.082
0.000	0.051	0.078	0.197	0.082
0.000	0.051	0.079	0.197	0.082
0.000	0.051	0.079	0.197	0.082
0.000	0.050	0.079	0.198	0.081
0.000	0.051	0.079	0.198	0.082
0.000	0.051	0.079	0.198	0.082
0.000	0.051	0.079	0.197	0.081
0.000	0.051	0.078	0.197	0.081
0.000	0.051	0.079	0.198	0.082


Standard Deviations	0.0003	0.0004	0.0005	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

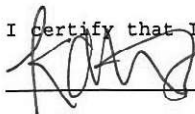

KATIE T SPEARIN

Signature and Printed Name

07/16/2026
Date



ANSI National Accreditation Board
ACCREDITED
ISO/IEC 17025
FORENSIC CALIBRATION
LABORATORY

Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-003945, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number: 80-003945
Owning Agency: LAKE LAND PD
Calibration Date: 07/16/2026
Calibration Time: 11:22

	UNCERTAINTY * $\pm$
0.050 g/ 210 L	0.004
0.080 g/ 210 L	0.004
0.200 g/ 210 L	0.008
0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.
Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).
The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/IEC 17025 standards.
This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Date

07/16/2026

FDLE/ATP Form 69 January 2026
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality

KATIE T SPEARIN,
Department Inspector