

Alcohol Testing Program - Instrument Processing Sheet

Agency: HILLSBOROUGH CO SO

Instrument Serial Number: 80-007482

Date In: 4/30/2026

DI Completion Date: 5/5/26

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By: <u>WKP</u> Date: <u>4/30/2026</u>		Quality Checks By: <u>IAS</u> Date: <u>5/1/2026</u>		Flow Column By: _____																																									
☒ Annual ☐ Dropped Off ☐ Registration ☒ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☐ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>215</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP102</u> 32 mm <u>.156</u> (.139-.169) 36 mm <u>.171</u> (.156-.190) 53 mm <u>.246</u> (.228-.278) 103 mm <u>.515</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>34419</u> Gauge: <u>1010</u> Instrument: <u>1012</u> ☒ Stability Checks		☐ Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)																																									
		<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot#/Exp</th></tr></thead><tbody><tr><td>0.050</td><td>MP5088</td><td>202603A 3/3/2028</td></tr><tr><td>0.080</td><td>MP5089</td><td>202603B 3/3/2028</td></tr><tr><td>0.200</td><td>MP5090</td><td>202603D 3/4/2028</td></tr><tr><td>0.080 DGS</td><td>23825080A1 N/A</td><td>23824080A1 10/5/2027</td></tr></tbody></table>		Simulator	Serial #	Lot#/Exp	0.050	MP5088	202603A 3/3/2028	0.080	MP5089	202603B 3/3/2028	0.200	MP5090	202603D 3/4/2028	0.080 DGS	23825080A1 N/A	23824080A1 10/5/2027	Maintenance By: _____ Date: _____ ☐ Battery Replacement ☐ Dry Gas Regulator Replacement ☐ Tank Sensor Tare ☐ Breath Tube Replacement ☐ Other:																										
Simulator	Serial #	Lot#/Exp																																											
0.050	MP5088	202603A 3/3/2028																																											
0.080	MP5089	202603B 3/3/2028																																											
0.200	MP5090	202603D 3/4/2028																																											
0.080 DGS	23825080A1 N/A	23824080A1 10/5/2027																																											
Optical Bench Adjustment By: _____		Department Inspection By: <u>KTS</u>																																											
Barometric Pressure Gauge: _____ ID#: _____		Barometric Pressure ID#: <u>28427</u>																																											
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.000</td><td></td><td>N/A</td><td>N/A</td></tr><tr><td>0.040</td><td></td><td></td><td></td></tr><tr><td>0.100</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.300</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge: <u>1010</u> Instrument: <u>1015</u> Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>9/25/2027</u> Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>9/22/2027</u>		<table border="1"><thead><tr><th>Simulator</th><th>Serial Number</th></tr></thead><tbody><tr><td>0.000</td><td>MP6289</td></tr><tr><td>Interferent</td><td>MP6290</td></tr><tr><td>0.050</td><td>MP6291</td></tr><tr><td>0.080</td><td>MP6292</td></tr><tr><td>0.200</td><td>MP6293</td></tr></tbody></table>		Simulator	Serial Number	0.000	MP6289	Interferent	MP6290	0.050	MP6291	0.080	MP6292	0.200	MP6293
Simulator	Serial #	Lot #	Expiration																																										
0.000		N/A	N/A																																										
0.040																																													
0.100																																													
0.200																																													
0.300																																													
0.080 DGS	N/A																																												
Simulator	Serial Number																																												
0.000	MP6289																																												
Interferent	MP6290																																												
0.050	MP6291																																												
0.080	MP6292																																												
0.200	MP6293																																												
☐ Post Optical Bench Adjustment Stability Checks		Attachments																																											
<table border="1"><thead><tr><th>Simulator</th><th>Serial #</th><th>Lot #</th><th>Expiration</th></tr></thead><tbody><tr><td>0.050</td><td></td><td></td><td></td></tr><tr><td>0.080</td><td></td><td></td><td></td></tr><tr><td>0.200</td><td></td><td></td><td></td></tr><tr><td>0.080 DGS</td><td>N/A</td><td></td><td></td></tr></tbody></table>		Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Optical Bench Adjustment		☐ Post-Stability Checks ☐ Flow Adjustment ☐ Form 40 ☐ Other:																					
Simulator	Serial #	Lot #	Expiration																																										
0.050																																													
0.080																																													
0.200																																													
0.080 DGS	N/A																																												
Gauge ID #: _____ Gauge: _____ Instrument: _____																																													
Notes/Suggested Service: Tech Review: Added comment to DI. KTS 5/18/26 Corrected DGS lot number. 5/19/25 KTS		☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																																											
		<table border="0"><tr><td>LeAndra Higginbotham</td><td>Digitally signed by LeAndra Higginbotham Date: 2026.05.20 09:42:00 -04'00'</td><td>Wen-Chi Pierson</td><td>Digitally signed by Wen-Chi Pierson Date: 2026.05.22 13:07:47 -04'00'</td></tr><tr><td colspan="2">Tech Review</td><td colspan="2">Admin Review</td></tr></table>				LeAndra Higginbotham	Digitally signed by LeAndra Higginbotham Date: 2026.05.20 09:42:00 -04'00'	Wen-Chi Pierson	Digitally signed by Wen-Chi Pierson Date: 2026.05.22 13:07:47 -04'00'	Tech Review		Admin Review																																	
LeAndra Higginbotham	Digitally signed by LeAndra Higginbotham Date: 2026.05.20 09:42:00 -04'00'	Wen-Chi Pierson	Digitally signed by Wen-Chi Pierson Date: 2026.05.22 13:07:47 -04'00'																																										
Tech Review		Admin Review																																											

Stability Checks

0.050 g/210L	0.080 g/210L	0.200 g/210L	DGS 0.080 g/210L																																																																																																																																																
0.047 to 0.053 g/210L ☒	0.077 to 0.083 g/210L ☒	0.194 to 0.206 g/210L ☒	0.077 to 0.083 g/210L ☒ ≤0.003 g/210L of Wet ☒																																																																																																																																																
HILLSBOROUGH CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007482 05/01/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:06</td></tr> <tr><td>Control Test</td><td>0.048</td><td>15:06</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:07</td></tr> <tr><td>Control Test</td><td>0.048</td><td>15:07</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:08</td></tr> <tr><td>Control Test</td><td>0.048</td><td>15:09</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:09</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0480</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	15:06	Control Test	0.048	15:06	Air Blank	0.000	15:07	Control Test	0.048	15:07	Air Blank	0.000	15:08	Control Test	0.048	15:09	Air Blank	0.000	15:09	Control Test Stats			Average	0.0480		Std Dev	0.0000		Rel Std Dev(%)	0.0000		HILLSBOROUGH CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007482 05/01/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>14:16</td></tr> <tr><td>Control Test</td><td>0.077</td><td>14:16</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:17</td></tr> <tr><td>Control Test</td><td>0.077</td><td>14:18</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:18</td></tr> <tr><td>Control Test</td><td>0.077</td><td>14:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>14:19</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0770</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	14:16	Control Test	0.077	14:16	Air Blank	0.000	14:17	Control Test	0.077	14:18	Air Blank	0.000	14:18	Control Test	0.077	14:19	Air Blank	0.000	14:19	Control Test Stats			Average	0.0770		Std Dev	0.0000		Rel Std Dev(%)	0.0000		HILLSBOROUGH CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007482 05/01/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:12</td></tr> <tr><td>Control Test</td><td>0.200</td><td>15:12</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:13</td></tr> <tr><td>Control Test</td><td>0.199</td><td>15:14</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:14</td></tr> <tr><td>Control Test</td><td>0.200</td><td>15:15</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:16</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1997</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2892</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	15:12	Control Test	0.200	15:12	Air Blank	0.000	15:13	Control Test	0.199	15:14	Air Blank	0.000	15:14	Control Test	0.200	15:15	Air Blank	0.000	15:16	Control Test Stats			Average	0.1997		Std Dev	0.0006		Rel Std Dev(%)	0.2892		HILLSBOROUGH CSO Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-007482 05/01/2026 Software: 8100.27 <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>15:18</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:19</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:19</td></tr> <tr><td>Control Test</td><td>0.079</td><td>15:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:20</td></tr> <tr><td>Control Test</td><td>0.080</td><td>15:20</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>15:21</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0793</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7277</td><td></td></tr> </tbody> </table> <i>Isabella Smykowski</i> Operator's Signature	Test	g/210L	Time	Air Blank	0.000	15:18	Control Test	0.079	15:19	Air Blank	0.000	15:19	Control Test	0.079	15:20	Air Blank	0.000	15:20	Control Test	0.080	15:20	Air Blank	0.000	15:21	Control Test Stats			Average	0.0793		Std Dev	0.0006		Rel Std Dev(%)	0.7277	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:06																																																																																																																																																	
Control Test	0.048	15:06																																																																																																																																																	
Air Blank	0.000	15:07																																																																																																																																																	
Control Test	0.048	15:07																																																																																																																																																	
Air Blank	0.000	15:08																																																																																																																																																	
Control Test	0.048	15:09																																																																																																																																																	
Air Blank	0.000	15:09																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0480																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	14:16																																																																																																																																																	
Control Test	0.077	14:16																																																																																																																																																	
Air Blank	0.000	14:17																																																																																																																																																	
Control Test	0.077	14:18																																																																																																																																																	
Air Blank	0.000	14:18																																																																																																																																																	
Control Test	0.077	14:19																																																																																																																																																	
Air Blank	0.000	14:19																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0770																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:12																																																																																																																																																	
Control Test	0.200	15:12																																																																																																																																																	
Air Blank	0.000	15:13																																																																																																																																																	
Control Test	0.199	15:14																																																																																																																																																	
Air Blank	0.000	15:14																																																																																																																																																	
Control Test	0.200	15:15																																																																																																																																																	
Air Blank	0.000	15:16																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1997																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2892																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	15:18																																																																																																																																																	
Control Test	0.079	15:19																																																																																																																																																	
Air Blank	0.000	15:19																																																																																																																																																	
Control Test	0.079	15:20																																																																																																																																																	
Air Blank	0.000	15:20																																																																																																																																																	
Control Test	0.080	15:20																																																																																																																																																	
Air Blank	0.000	15:21																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0793																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7277																																																																																																																																																		

Florida Department of Law Enforcement

Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: HILLSBOROUGH CSO
Time of Inspection: 10:16

Date of Inspection: 05/05/2026

Serial Number: 80-007482
Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202603D Exp: 03/04/2028	0.08 g/210L Dry Gas Std Test (g/210L) 23225080A1 Lot#:202603D Exp: 03/04/2028
0.000	0.049	0.078	0.201	0.080
0.000	0.049	0.078	0.201	0.080
0.000	0.049	0.077	0.201	0.080
0.000	0.049	0.078	0.201	0.080
0.000	0.049	0.078	0.201	0.080
0.000	0.049	0.078	0.201	0.080
0.000	0.049	0.078	0.201	0.080
0.000	0.049	0.078	0.201	0.079
0.000	0.049	0.078	0.201	0.079
0.000	0.049	0.078	0.201	0.080

KTS
5/5/26

Standard Deviations	0.0000	0.0003	0.0000	0.0004
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0001 Number of Simulators Used: 5


Remarks:

A F / M A: No Sample Provided.

Test timed out before sample was provided. Repeated test and provided sample in allowable time. KTS 5/18/26

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

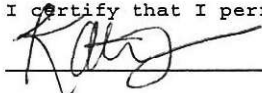

KATIE T SPEARIN

Signature and Printed Name

05/05/2026
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007482, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007482</u>	UNCERTAINTY* $\pm$	
Owning Agency:	<u>HILLSBOROUGH CSO</u>	0.050 g/ 210 L	0.004
Calibration Date:	<u>05/05/2026</u>	0.080 g/ 210 L	0.004
Calibration Time:	<u>10:16</u>	0.200 g/ 210 L	0.008
		0.080 g/ 210 L Dry Gas Control	0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.


Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

05/05/2026

Date

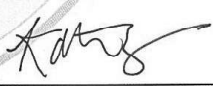
KATIE T SPEARIN,
Department Inspector