

Alcohol Testing Program - Instrument Processing Sheet

Agency: GROVELAND PD

Instrument Serial Number: 80-007252

Date In: 8/4/2026

DI Completion Date: 8/10/26

☒ Ship ☐ P/U ☐ H/D ☐ CMI ☐ EE

Intake By: <u>WKP</u> Date: <u>8/5/2026</u>		Quality Checks By: <u>KTS</u> Date: <u>8/6/2026</u>		Flow Adjustment By: _____			
☒ Annual ☐ Dropped Off ☐ Registration ☐ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☐ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>168</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP103</u> 32 mm <u>0.160</u> (.139-.169) 36 mm <u>0.175</u> (.156-.190) 53 mm <u>0.246</u> (.228-.278) 103 mm <u>0.488</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1015</u> Instrument: <u>1014</u> ☒ Stability Checks		Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)			
						Maintenance By: _____ Date: _____	
						☐ Battery Replacement	
						☐ Dry Gas Regulator Replacement	
						☐ Tank Sensor Tare	
☐ Breath Tube Replacement							
☐ Other:							
Optical Bench Adjustment By: _____		Department Inspection By: <u>KTS</u>					
Barometric Pressure Gauge: _____ ID#: _____		Barometric Pressure ID#: <u>28427</u>					
Simulator	Serial #	Lot #	Expiration	Gauge: <u>1014</u>	Instrument: <u>1014</u>		
0.000		N/A	N/A	Mouth Alcohol Solution Lot #: <u>2025-D</u>	Exp: <u>9/25/27</u>		
0.040				Acetone Stock Solution Lot #: <u>2025-B</u>	Exp: <u>9/22/27</u>		
0.100				Simulator	Serial Number		
0.200				0.000	MP5086		
0.300				Interferent	MP5087		
0.080 DGS	N/A			0.050	MP5088		
☐ Post Optical Bench Adjustment Stability Checks				0.080	MP5089		
Simulator	Serial #	Lot #	Expiration	0.200	MP5090		
0.050				Attachments			
0.080				☒ Form 41			
0.200				☒ Stability Checks			
0.080 DGS	N/A			☒ Calibration Certificate			
Gauge ID #: _____				☐ Optical Bench Adjustment			
Gauge: _____ Instrument: _____				☐ Post-Stability Checks			
Notes/Suggested Service:				☐ Flow Adjustment			
				☒ Form 40			
				☒ Other: Calibration Certificate Calculations			
				☒ Instrument Complies with Chapter 11D-8, FAC			
				☐ Instrument Does Not Comply with Chapter 11D-8, FAC			
				☒ Return to/Place into Evidentiary Use			
				☐ Remain Out of Evidentiary Use			
				☒ Conduct an Agency Inspection Before Evidentiary Use			
Wen-Chi Pierson		Digitally signed by Wen-Chi Pierson Date: 2026.08.13 14:39:37 -04'00'		S. LeAndra Higginbotham			
				Digitally signed by S. LeAndra Higginbotham Date: 2026.08.14 10:15:51 -04'00'			
Tech Review				Admin Review			

Florida Department of Law Enforcement Alcohol Testing Program

AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: GROVELAND PD

Time of Inspection: 06:26

Date of Inspection: 08/06/2026

Serial Number: 80-007252

Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: Exp:	0.08g/210L Test (g/210L) Lot#: Exp:	0.20g/210L Test (g/210L) Lot#: Exp:	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: Exp:


Number of Simulators Used: _____

Remarks:

BYPASSED TO OPERATE- COMPLIANCE NOT DETERMINED

The above instrument complies (X) does not comply () with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

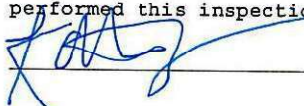


KATIE T SPEARIN
Signature and Printed Name

08/06/2026
Date

DG S

GROVELAND PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-007252
08/06/2026
Software: 8100.27


Test	g/210L	Time
Air Blank	0.000	07:03
Control Test	0.080	07:04
Air Blank	0.000	07:05
Control Test	0.079	07:05
Air Blank	0.000	07:06
Control Test	0.080	07:07
Air Blank	0.000	07:07
Control Test Stats		
Average	0.0797	
Std Dev	0.0006	
Rel Std Dev(%)	0.7247	

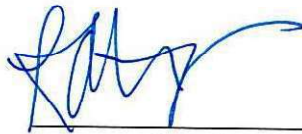

Operator's Signature

GROVELAND PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-007252
08/06/2026
Software: 8100.27


Test	g/210L	Time
Air Blank	0.000	06:58
Control Test	0.050	06:59
Air Blank	0.000	07:00
Control Test	0.051	07:00
Air Blank	0.000	07:01
Control Test	0.050	07:01
Air Blank	0.000	07:02
Control Test Stats		
Average	0.0503	
Std Dev	0.0006	
Rel Std Dev(%)	1.1471	

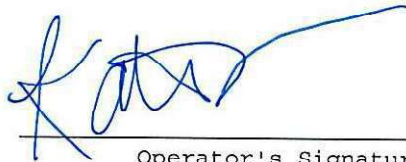

Operator's Signature

GROVELAND PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-007252
08/06/2026
Software: 8100.27


Test	g/210L	Time
Air Blank	0.000	06:52
Control Test	0.079	06:53
Air Blank	0.000	06:53
Control Test	0.079	06:53
Air Blank	0.000	06:54
Control Test	0.079	06:54
Air Blank	0.000	06:55
Control Test Stats		
Average	0.0790	
Std Dev	0.0000	
Rel Std Dev(%)	0.0000	

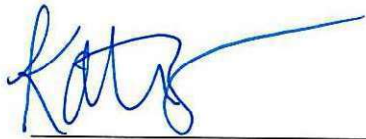

Operator's Signature

GROVELAND PD
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-007252
08/06/2026
Software: 8100.27

Test	g/210L	Time
Air Blank	0.000	07:09
Control Test	0.202	07:09
Air Blank	0.000	07:10
Control Test	0.202	07:11
Air Blank	0.000	07:11
Control Test	0.201	07:12
Air Blank	0.000	07:13
Control Test Stats		
Average	0.2017	
Std Dev	0.0006	
Rel Std Dev(%)	0.2863	

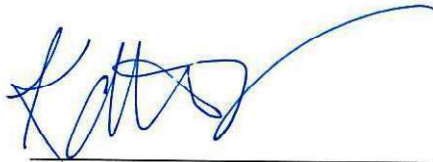

Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: GROVELAND PD

Time of Inspection: 09:00

Date of Inspection: 08/10/2026

Serial Number: 80-007252

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202603D Exp: 03/04/2028	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:23825080A1 Exp: 10/05/2027
0.000	0.049	0.080	0.202	0.080
0.000	0.050	0.080	0.202	0.079
0.000	0.050	0.080	0.202	0.080
0.000	0.050	0.080	0.202	0.080
0.000	0.050	0.080	0.202	0.079
0.000	0.050	0.081	0.203	0.079
0.000	0.051	0.081	0.202	0.079
0.000	0.050	0.081	0.202	0.080
0.000	0.051	0.081	0.202	0.080
0.000	0.050	0.081	0.202	0.079

Standard Deviations	0.0005	0.0005	0.0003	0.0005
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0004 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

KATIE T SPEARIN

Signature and Printed Name

08/10/2026
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-007252, manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-007252</u>	Average of y ± Uncertainty*
Owning Agency:	<u>GROVELAND PD</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>08/10/2026</u>	0.081 g/ 210 L 0.004
Calibration Time:	<u>09:00</u>	0.202 g/ 210 L 0.008
		0.080 g/ 210 L Dry Gas Control 0.005

All results are reported in g/210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.
*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence (k=3).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Digitally signed by Kaitlyn
Spearin
Date: 2026.08.11 12:54:58
-04'00'

Kaitlyn Spearin

08/10/2026

Date

KATIE T SPEARIN,

Department Inspector

FDLE/ATP Form 69 August 2026
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality