

Alcohol Testing Program - Instrument Processing Sheet

Agency: FL HIGHWAY PATROL Instrument Serial Number: 80-006765  
 Date In: 2/9/2026 DI Completion Date: N/A SLH  Ship  P/U  H/D  CMI  EE

Intake By: <u>WKP</u> Date: <u>2/9/2026</u> <input checked="" type="checkbox"/> Annual <input type="checkbox"/> Dropped Off <input type="checkbox"/> Registration <input type="checkbox"/> Return from CMI / EE <input type="checkbox"/> Training Instrument Visual Inspection <input checked="" type="checkbox"/> Case <input checked="" type="checkbox"/> Handle <input checked="" type="checkbox"/> Keyboard <input checked="" type="checkbox"/> Dry Gas Shelf <input checked="" type="checkbox"/> Feet <input checked="" type="checkbox"/> Breath Tube <input checked="" type="checkbox"/> Ports <input checked="" type="checkbox"/> Screws Tight Other Equipment/Accessories <input type="checkbox"/> Power Cord <input type="checkbox"/> Printer Cable <input checked="" type="checkbox"/> Static Bag <input type="checkbox"/> 12V DC Cable Notes: Portion of a possible breath test with an analytical stability failed diagnostic was taped to the instrument. Contacted agency for additional information. SLH 2/9/2026	Quality Checks By: <u>SLH</u> Date: <u>2/11/2026</u> <input checked="" type="checkbox"/> Breath Tube Screen <input checked="" type="checkbox"/> Replace External O-Rings <input checked="" type="checkbox"/> Instrument Set Up Verified <input checked="" type="checkbox"/> R-Value: <u>220</u> <input checked="" type="checkbox"/> Flow Verification (L/s) Flow Column #: <u>ATP 102</u> 32 mm <u>0.144</u> (.139-.169) 36 mm <u>0.160</u> (.156-.190) 53 mm <u>0.234</u> (.228-.278) 103 mm <u>0.488</u> (.447-.547) <input checked="" type="checkbox"/> Barometric Pressure Check Gauge ID #: <u>28421</u> Gauge: <u>1019</u> Instrument: <u>1017</u> <input checked="" type="checkbox"/> Stability Checks	Flow Adjustment By: _____ Flow Column #: _____ <input type="checkbox"/> 5L/min – 17mm <input type="checkbox"/> 15L/min – 53mm <input type="checkbox"/> 30L/min – 103mm <input type="checkbox"/> R-Value: _____ <input type="checkbox"/> Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)																				
<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <td style="width:25%;"></td> <td style="width:25%;"></td> <td style="width:25%;"></td> <td style="width:25%;"></td> </tr> <tr> <td>0.050</td> <td>MP6291</td> <td>202406K</td> <td>06/19/2026</td> </tr> <tr> <td>0.080</td> <td>MP6292</td> <td>202406L</td> <td>06/19/2026</td> </tr> <tr> <td>0.200</td> <td>MP6293</td> <td>202406N</td> <td>06/20/2026</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>AG510701</td> <td>04/17/2027</td> </tr> </table>							0.050	MP6291	202406K	06/19/2026	0.080	MP6292	202406L	06/19/2026	0.200	MP6293	202406N	06/20/2026	0.080 DGS	N/A	AG510701	04/17/2027
0.050	MP6291	202406K	06/19/2026																			
0.080	MP6292	202406L	06/19/2026																			
0.200	MP6293	202406N	06/20/2026																			
0.080 DGS	N/A	AG510701	04/17/2027																			
Maintenance By: <u>SLH</u> Date: <u>2/12/2026</u> <input type="checkbox"/> Battery Replacement <input type="checkbox"/> Dry Gas Regulator Replacement <input type="checkbox"/> Tank Sensor Tare <input type="checkbox"/> Breath Tube Replacement <input checked="" type="checkbox"/> Other: <b>Added internal printer paper</b>																						

<b>Optical Bench Adjustment</b> By: _____	<b>Department Inspection</b> By: _____																																								
Barometric Pressure Gauge: _____ ID#: _____	Barometric Pressure ID#: _____																																								
<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> <tr> <td>0.000</td> <td></td> <td>N/A</td> <td>N/A</td> </tr> <tr> <td>0.040</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.100</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.300</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>	Simulator	Serial #	Lot #	Expiration	0.000		N/A	N/A	0.040				0.100				0.200				0.300				0.080 DGS	N/A			Gauge: _____ Instrument: _____ Mouth Alcohol Solution Lot #: _____ Exp: _____ Acetone Stock Solution Lot #: _____ Exp: _____ <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial Number</th> </tr> <tr> <td>0.000</td> <td></td> </tr> <tr> <td>Interferent</td> <td></td> </tr> <tr> <td>0.050</td> <td></td> </tr> <tr> <td>0.080</td> <td></td> </tr> <tr> <td>0.200</td> <td></td> </tr> </table>	Simulator	Serial Number	0.000		Interferent		0.050		0.080		0.200	
Simulator	Serial #	Lot #	Expiration																																						
0.000		N/A	N/A																																						
0.040																																									
0.100																																									
0.200																																									
0.300																																									
0.080 DGS	N/A																																								
Simulator	Serial Number																																								
0.000																																									
Interferent																																									
0.050																																									
0.080																																									
0.200																																									
<input type="checkbox"/> Post Optical Bench Adjustment Stability Checks																																									
<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot #</th> <th>Expiration</th> </tr> <tr> <td>0.050</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.200</td> <td></td> <td></td> <td></td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td></td> <td></td> </tr> </table>	Simulator	Serial #	Lot #	Expiration	0.050				0.080				0.200				0.080 DGS	N/A			<b>Attachments</b> <input type="checkbox"/> Form 41 <input checked="" type="checkbox"/> Stability Checks <input type="checkbox"/> Calibration Certificate <input type="checkbox"/> Optical Bench Adjustment <input type="checkbox"/> Post-Stability Checks <input type="checkbox"/> Flow Adjustment <input checked="" type="checkbox"/> Form 40 <input checked="" type="checkbox"/> Other: Portion of a test submitted by agency taped to the instrument, Diagnostic tests performed by DI (x2), and RMA																				
Simulator	Serial #	Lot #	Expiration																																						
0.050																																									
0.080																																									
0.200																																									
0.080 DGS	N/A																																								
Gauge ID #: _____ Gauge: _____ Instrument: _____																																									

Notes/Suggested Service: RMA = return material authorization. SLH 2/17/2026 Due to failed diagnostic submitted by agency and discussion with CMI on 2/12/26, will send for repair. SLH 2/17/26 Stability check performed 2/12/26. SLH 2/12/2026	<input type="checkbox"/> Instrument Complies with Chapter 11D-8, FAC <input type="checkbox"/> Instrument Does Not Comply with Chapter 11D-8, FAC <input type="checkbox"/> Return to/Place into Evidentiary Use <input checked="" type="checkbox"/> Remain Out of Evidentiary Use <input type="checkbox"/> Conduct an Agency Inspection Before Evidentiary Use <div style="text-align: right;">                 Digitally signed by Shayla Platt                  Date: 2026.02.20 23:01:03 -05'00'                  Digitally signed by Shayla Platt                  Date: 2026.02.20 23:00:37 -05'00'             </div>
Shayla Platt Tech Review	Shayla Platt Admin Review

80-006765- this paper portion of a possible  
breath test was taped to the instrument when  
submitted. Although illegible in the scan read  
Analytical Stability Test diagnostic failure.  
SLH 2/9/2026

Test	g/210L	Time
Diagnostics Check	OK	00:39
Air Blank	0.000	00:39
Control Test	0.083	00:39
Air Blank	0.000	00:40
Subject Sample #1	0.228	00:41
Air Blank	0.000	00:41
Air Blank	0.000	00:43
Subject Sample #2	0.212	00:43
Air Blank	0.000	00:44
Control Test	0.077	00:44
Air Blank	0.000	00:45
Diagnostics Check	Fail	00:45
Air Blank	0.000	00:46

Analytical Stability Test

FL HIGHWAY PATROL  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-006765  
02/12/2026  
Software: 8100.27

-----  
: DIAGNOSTICS :  
-----

Voltage/Current Test	OK
RAM Test	OK
EEPROM Checksum Test	OK
Real Time Clock Test	OK
DSP Test	OK
Analytical Stability Test	OK
Internal Printer Test	OK
Modem Test	OK
Temperature Regulation Test	OK

FL HIGHWAY PATROL  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-006765  
02/12/2026  
Software: 8100.27

-----  
: DIAGNOSTICS :  
-----

Voltage/Current Test	OK
RAM Test	OK
EEPROM Checksum Test	OK
Real Time Clock Test	OK
DSP Test	OK
Analytical Stability Test	OK
Internal Printer Test	OK
Modem Test	OK
Temperature Regulation Test	OK

Analytical Stability Test  
passed diagnostics (x2).

SUA

2/12/2026

# Florida Department of Law Enforcement Alcohol Testing Program

## AGENCY INSPECTION REPORT - INTOXILYZER 8000

Agency: FL HIGHWAY PATROL  
Time of Inspection: 10:47

Date of Inspection: 02/11/2026

Serial Number: 80-006765  
Software: 8100.27

Check or Test	YES	NO
Date and/or Time Adjusted		No
Diagnostic Check (Pre-Inspection): OK		No
Alcohol Free Subject Test: 0.000		No
Mouth Alcohol Test: Slope Not Met		No
Interferent Detect Test: Interferent Detect		No
Diagnostic Check (Post-Inspection): OK		No


Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08g/210L Test (g/210L) Lot#: _____ Exp: _____	0.20g/210L Test (g/210L) Lot#: _____ Exp: _____	0.08 g/210L Dry Gas Std Test (g/210L) Lot#: _____ Exp: _____

Number of Simulators Used: \_\_\_\_\_

Remarks:  
BYPASS AI FOR OPERATION, COMPLIANCE UNDETERMINED

The above instrument complies (  ) does not comply (  ) with Chapter 11D-8, FAC.

I certify that I hold a valid Florida Department of Law Enforcement Agency Inspector Permit and that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

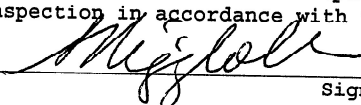

LEANDRA HIGGINBOTHAM  
Signature and Printed Name

02/11/2026  
Date

# Stability Checks

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L SO.003 g/210L of Wet																																																																																																																																																
<p>FL HIGHWAY PATROL Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006765 02/12/2026 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>11:57</td></tr> <tr><td>Control Test</td><td>0.048</td><td>11:58</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>11:58</td></tr> <tr><td>Control Test</td><td>0.049</td><td>11:59</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:00</td></tr> <tr><td>Control Test</td><td>0.049</td><td>12:00</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>12:01</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.1487</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.1853</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	11:57	Control Test	0.048	11:58	Air Blank	0.000	11:58	Control Test	0.049	11:59	Air Blank	0.000	12:00	Control Test	0.049	12:00	Air Blank	0.000	12:01	Control Test Stats			Average	0.1487		Std Dev	0.0006		Rel Std Dev(%)	0.1853		<p>FL HIGHWAY PATROL Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006765 02/12/2026 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:51</td></tr> <tr><td>Control Test</td><td>0.079</td><td>10:51</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:52</td></tr> <tr><td>Control Test</td><td>0.080</td><td>10:53</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:53</td></tr> <tr><td>Control Test</td><td>0.080</td><td>10:54</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:54</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0797</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.7247</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:51	Control Test	0.079	10:51	Air Blank	0.000	10:52	Control Test	0.080	10:53	Air Blank	0.000	10:53	Control Test	0.080	10:54	Air Blank	0.000	10:54	Control Test Stats			Average	0.0797		Std Dev	0.0006		Rel Std Dev(%)	0.7247		<p>FL HIGHWAY PATROL Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006765 02/12/2026 Software: 8110.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>10:33</td></tr> <tr><td>Control Test</td><td>0.203</td><td>10:34</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:35</td></tr> <tr><td>Control Test</td><td>0.202</td><td>10:35</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:36</td></tr> <tr><td>Control Test</td><td>0.202</td><td>10:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>10:37</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.2123</td><td></td></tr> <tr><td>Std Dev</td><td>0.0006</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.2853</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	10:33	Control Test	0.203	10:34	Air Blank	0.000	10:35	Control Test	0.202	10:35	Air Blank	0.000	10:36	Control Test	0.202	10:37	Air Blank	0.000	10:37	Control Test Stats			Average	0.2123		Std Dev	0.0006		Rel Std Dev(%)	0.2853		<p>FL HIGHWAY PATROL Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-006765 02/12/2026 Software: 8100.27</p> <table border="1"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.000</td><td>09:36</td></tr> <tr><td>Control Test</td><td>0.080</td><td>09:36</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>09:37</td></tr> <tr><td>Control Test</td><td>0.080</td><td>09:37</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>09:38</td></tr> <tr><td>Control Test</td><td>0.080</td><td>09:38</td></tr> <tr><td>Air Blank</td><td>0.000</td><td>09:38</td></tr> <tr><td>Control Test Stats</td><td></td><td></td></tr> <tr><td>Average</td><td>0.0800</td><td></td></tr> <tr><td>Std Dev</td><td>0.0000</td><td></td></tr> <tr><td>Rel Std Dev(%)</td><td>0.0000</td><td></td></tr> </tbody> </table>	Test	g/210L	Time	Air Blank	0.000	09:36	Control Test	0.080	09:36	Air Blank	0.000	09:37	Control Test	0.080	09:37	Air Blank	0.000	09:38	Control Test	0.080	09:38	Air Blank	0.000	09:38	Control Test Stats			Average	0.0800		Std Dev	0.0000		Rel Std Dev(%)	0.0000	
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	11:57																																																																																																																																																	
Control Test	0.048	11:58																																																																																																																																																	
Air Blank	0.000	11:58																																																																																																																																																	
Control Test	0.049	11:59																																																																																																																																																	
Air Blank	0.000	12:00																																																																																																																																																	
Control Test	0.049	12:00																																																																																																																																																	
Air Blank	0.000	12:01																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.1487																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.1853																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:51																																																																																																																																																	
Control Test	0.079	10:51																																																																																																																																																	
Air Blank	0.000	10:52																																																																																																																																																	
Control Test	0.080	10:53																																																																																																																																																	
Air Blank	0.000	10:53																																																																																																																																																	
Control Test	0.080	10:54																																																																																																																																																	
Air Blank	0.000	10:54																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0797																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.7247																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	10:33																																																																																																																																																	
Control Test	0.203	10:34																																																																																																																																																	
Air Blank	0.000	10:35																																																																																																																																																	
Control Test	0.202	10:35																																																																																																																																																	
Air Blank	0.000	10:36																																																																																																																																																	
Control Test	0.202	10:37																																																																																																																																																	
Air Blank	0.000	10:37																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.2123																																																																																																																																																		
Std Dev	0.0006																																																																																																																																																		
Rel Std Dev(%)	0.2853																																																																																																																																																		
Test	g/210L	Time																																																																																																																																																	
Air Blank	0.000	09:36																																																																																																																																																	
Control Test	0.080	09:36																																																																																																																																																	
Air Blank	0.000	09:37																																																																																																																																																	
Control Test	0.080	09:37																																																																																																																																																	
Air Blank	0.000	09:38																																																																																																																																																	
Control Test	0.080	09:38																																																																																																																																																	
Air Blank	0.000	09:38																																																																																																																																																	
Control Test Stats																																																																																																																																																			
Average	0.0800																																																																																																																																																		
Std Dev	0.0000																																																																																																																																																		
Rel Std Dev(%)	0.0000																																																																																																																																																		
<p>Operator: <i>M. J. J. J.</i> Signature</p>	<p>Operator: <i>M. J. J. J.</i> Signature</p>	<p>Operator: <i>M. J. J. J.</i> Signature</p>	<p>Operator: <i>M. J. J. J.</i> Signature</p>																																																																																																																																																

Return Material Authorization

**Ship to:**

CMI, Inc.

Enforcement Electronics

Shipment to repair facility authorized by: Joseph Farley on 2/9/2026

Items Returned:  Instrument  Supplies  Other Describe: \_\_\_\_\_

Instrument Model: Intoxilyzer 8000 Serial Number: 80-006765

Bill To Address:  Sergeant Joseph Farley Criminal Interdiction Unit Florida Highway Patrol 7322 Normandy Boulevard Jacksonville, Florida 32205	Ship to Address:  FDLE Offsite Mail Facility c/o FDLE HQ 813 B Lake Bradford Rd Tallahassee, FL 32304
------------------------------------------------------------------------------------------------------------------------------------------------------------------	----------------------------------------------------------------------------------------------------------------------

Reason for Return:  
During the DVM monitoring of the span noise, there was consistent large jumping of the counts on the 3um detector greater than 20 counts over several hours of observation.

I require an estimate **BEFORE** any repairs will be authorized and/or conducted

Please contact: Sgt Joseph Farley Phone #: 850-672-8222

Email: josephfarley@flhsmv.gov

ATP Contact Name: LeAndra Higginbotham ATP Email: leandrahigginbotham@fdle.state.fl.us